

# HPLC Application

ID No.: 15538

## Patulin analysis using Synergi Max-RP

**Column:** Synergi<sup>TM</sup> 4  $\mu$ m Max-RP 80 Å, LC Column 150 x 3 mm, Ea

**Dimensions:** 150 x 3 mm ID

**Order No:** 00F-4337-Y0

**Elution Type:** Isocratic

**Eluent A:** Water

**Eluent B:** Acetonitrile

| Gradient | Step No. | Time (min) | Pct A | Pct B |
|----------|----------|------------|-------|-------|
| Profile: | 1        | 0          | 90    | 10    |

**Flow Rate:** 0.5 mL/min

**Col. Temp.:** ambient

**Detection:** Refractive Index @ 0.0000000000 (ambient)


**Analyst Note:** SecurityGuard<sup>TM</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 2.0mm, 10/Pk Part No.: AJ0-6073

- Holder Part No.: KJ0-4282



15538

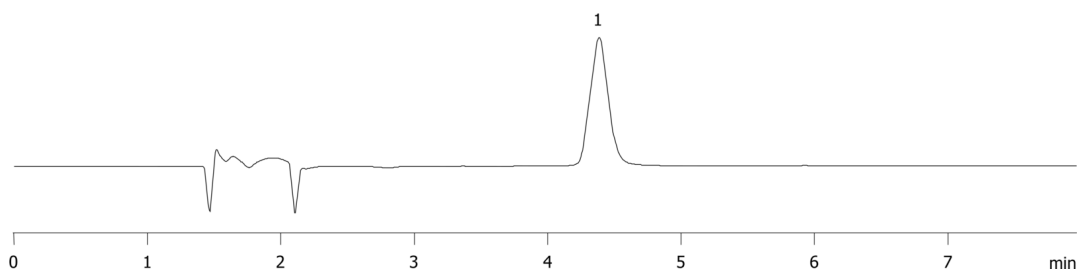

### ANALYTES:

1 Patulin

