


Alternative Polynuclear Aromatic Hydrocarbons (PAH) on ZB-5MS**Column:** Zebtron™ ZB-5MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 0.18 µm**Order No:** 7FD-G010-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 15.0000000000:1 0 µL @ 285°C**Detection:** Refractive Index (180°C)

15561

**Products used in this application:****ANALYTES:**

- 1 Benzo[j]fluoranthene
- 2 3-Methylcholanthrene
- 3 Dibenz[a,h]acridine
- 4 Dibenz[a,,j]acridine
- 5 7H-Dibenzo[c,g]carbazole
- 6 Dibenzo[a,e]pyrene
- 7 Dibenzo[a,i]pyrene
- 8 Dibenzo[a,j]pyrene

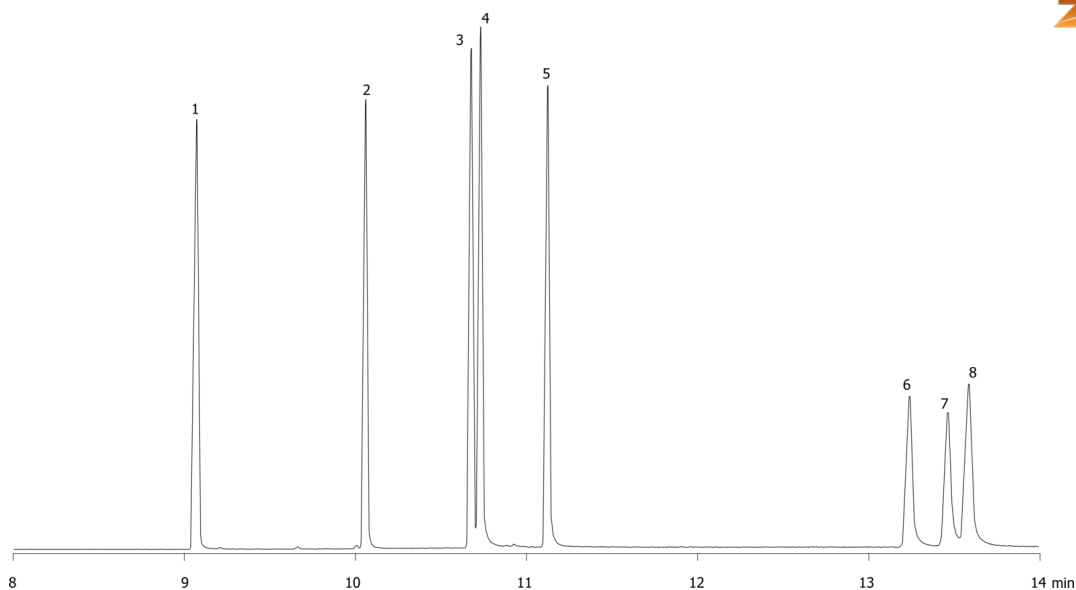