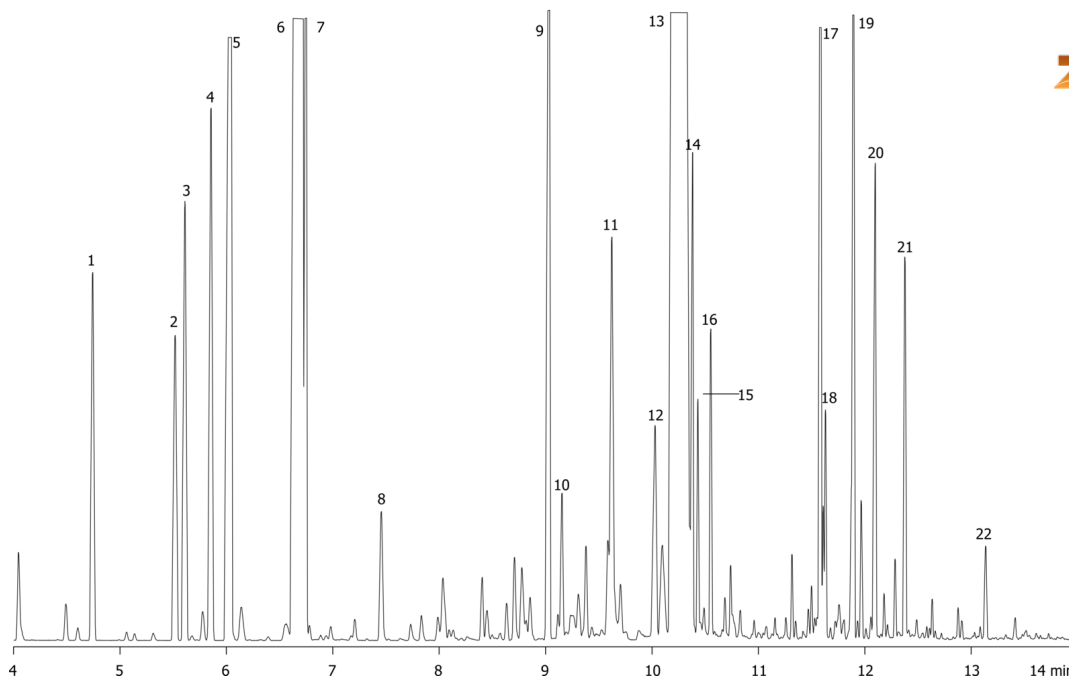


Spearmint Oil on ZB5ms 20x0.18x0.18**Column:** Zebron[™] ZB-5MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 0.18 µm**Order No:** 7FD-G010-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 100.000000000:1 0 µL @ 180°C**Detection:** Refractive Index (180°C)

15565

**Products used in this application:**

Spearmint Oil on ZB5ms 20x0.18x0.18

ANALYTES:

- 1 alpha-Pinene
- 2 Sabinene
- 3 beta-Pinene
- 4 beta-Myrcene
- 5 3-Octanal
- 6 Limonene
- 7 Eucalyptol
- 8 beta-Terpineol
- 9 Menthone
- 10 Menthone
- 11 Dihydrocarvone
- 12 Carveol
- 13 Carvone
- 14 Piperitone
- 15 trans-Carvone oxide
- 16 cis-Carvone Oxide
- 17 Bourbonene
- 18 Jasmone
- 19 Caryophyllene
- 20 Farnesene
- 21 Germacrene
- 22 Caryophyllene oxide

