


**Amphetamines (PFPA Derivatization) on ZB1ms 20x0.18x0.18****Column:** Zebtron<sup>™</sup> ZB-1MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 0.18 µm**Order No:** 7FD-G011-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 290°C**Detection:** Refractive Index (240°C)

15574

**Products used in this application:****ANALYTES:**

- 1 Amphetamine (PFP)
- 2 Methamphetamine-PFP

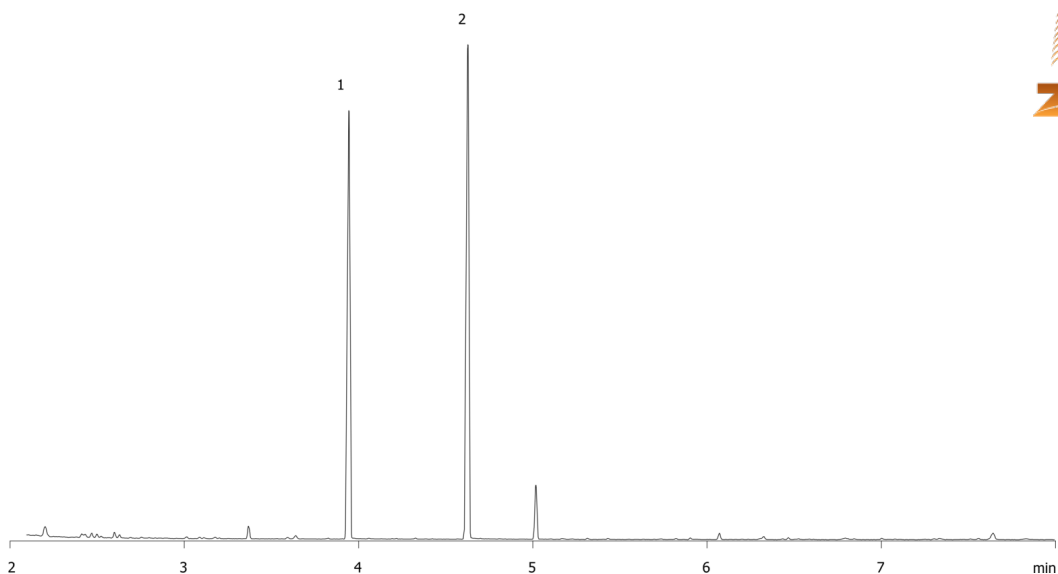