

Estradiol - USP Method (Synergi 4µm Fusion-RP)

Column: Synergi™ 4 µm Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4424-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	45	55

Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: Refractive Index @ 0.0000000000 (25 °C)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282



15663

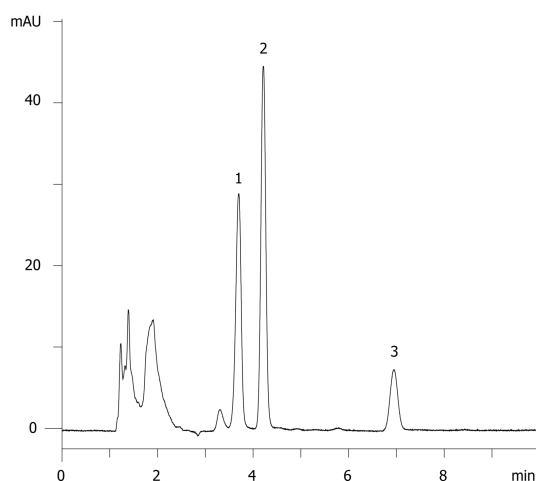

ANALYTES:

1 Ethyl paraben

Retention Time: 3.697 min

2 Estradiol

Retention Time: 4.213 min

3 Estrone

Retention Time: 6.946 min

