

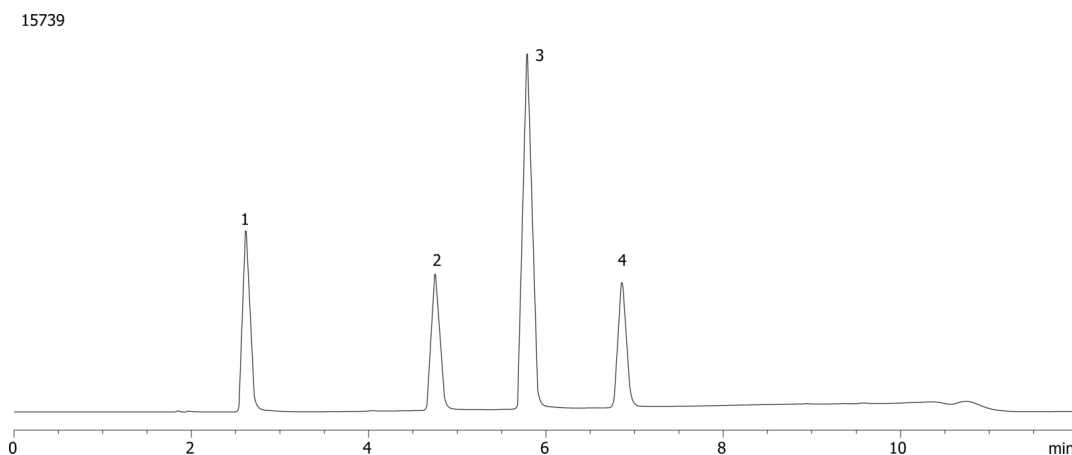
Catecholamines metabolites using Synergi Fusion-RP**Column:** SynergiTM 4 μ m Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4424-E0**Elution Type:** Gradient**Eluent A:** 20mM Ammonium formate, pH 3.5**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	10	50	50

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Vanillylmandelic acid
- 2 3,4-Dihydroxyphenylacetic acid (DOPAC)
- 3 5-Hydroxyindoleacetic acid
- 4 Homovanillic acid

