


Massachusetts EPH - Aromatic Fraction on ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** Splitless 0.000000000:1 0 μ L @ 285°C**Detection:** Refractive Index (315°C)**Analyst Note:** Analytes are 50ppm and IS and Surrogates are 40ppm for fractionation.**Products used in this application:**

15788

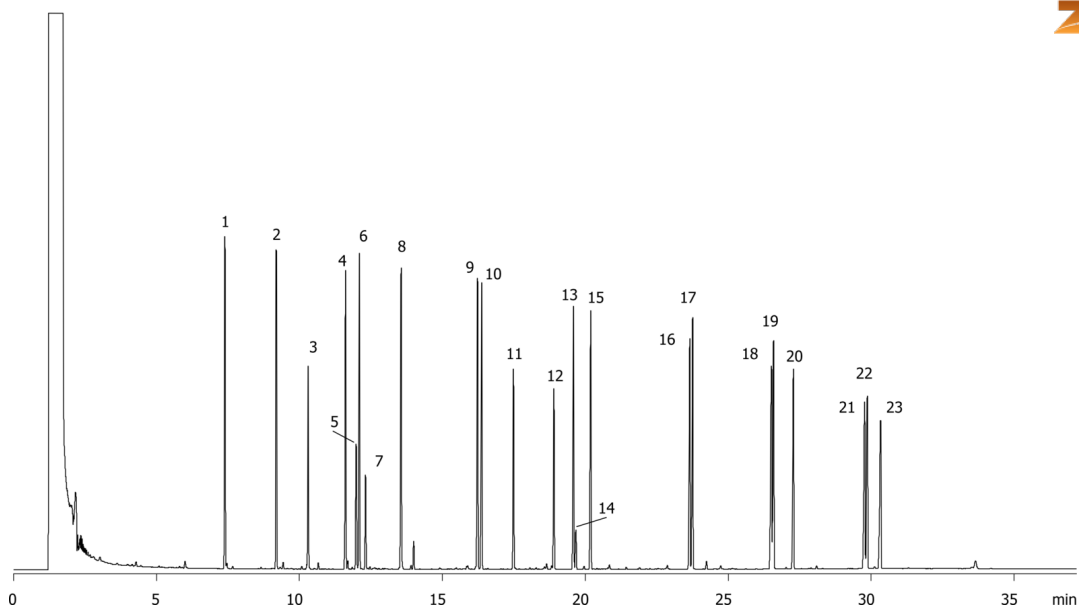

Massachusetts EPH - Aromatic Fraction on ZB-5MS**ANALYTES:**

- 1** Naphthalene
- 2** 2-Methylnaphthalene
- 3** 2-Fluorobiphenyl
- 4** Acenaphthalene
- 5** 2-Bromonaphthalene
- 6** Acenaphthene
- 7** phthalate
- 8** Fluorene
- 9** Phenanthrene
- 10** Anthracene
- 11** o-Terphenyl
- 12** 5-alpha-Androstane
- 13** Fluoranthene
- 14** 1-Chloro-octadecane
- 15** Pyrene
- 16** Benz[a]anthracene
- 17** Chrysene
- 18** Benzo[b]fluoranthene
- 19** Benzo[k]fluoranthene
- 20** Benzo[a]pyrene
- 21** Indeno[1,2,3-cd]pyrene
- 22** Dibenz[a,h]anthracene
- 23** Benzo[g,h,i]perylene

