


USP Glycerin on ZB-624**Column:** ZebtronTM ZB-624, GC Cap. Column 30 m x 0.53 mm x 3.00 μ m, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.53 mm x 3 μ m**Order No:** 7HK-G005-36**Oven Profile:****Carrier Gas:** Constant Flow Helium, 4.5 mL/min**Injection:** Split 10.0000000000:1 0 μ L @ 220°C**Detection:** Refractive Index (250°C)**Analyst Note:** Analytes are 0.5 mg/mL
Resolution (0.5mg/mL) 8.7
%RSD Diethylene Glycol (0.05mg/mL): 7.4%
%RSD Glycerine (0.05mg/mL): 6.9%**Products used in this application:****ANALYTES:**

- 1 Diethylene glycol
- 2 Glycerin

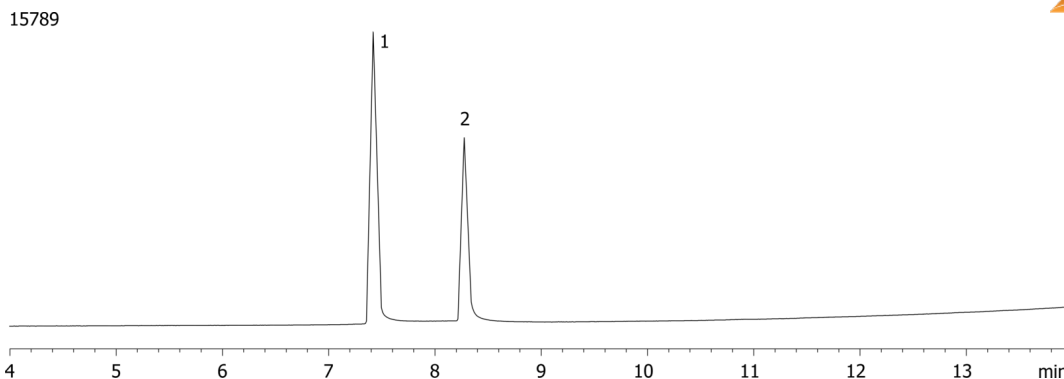