

## siRNA Cartridge Purification - Wash #1 Fraction

**Column:** Clarity® 3 µm Oligo-RP, LC Column 50 x 4.6 mm, Ea

**Dimensions:** 50 x 4.6 mm ID

**Order No:** 00B-4441-E0

**Elution Type:** Gradient

**Eluent A:** 50mM TEAA, pH 7.5 / 5% ACN

**Eluent B:** Methanol

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 100   | 0     |
|                   | 2        | 20         | 40    | 60    |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** Refractive Index @ 0.0000000000 (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Clarity Oligo-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-8135

- Holder Part No.: KJ0-4282

15802



Products used in this application:

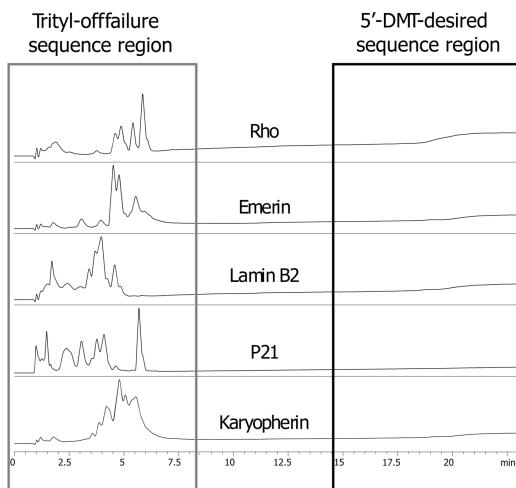

### ANALYTES:

- 1 Intital wash samples

