

## Drug Screen on Zebron ZB-1MS 20m x 0.18mm x 0.18um

**Column:** Zebron™ ZB-1MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea

**Phase:** 100% Dimethylpolysiloxane

**Dimensions:** 20 meters x 0.18 mm x 0.18 µm

**Order No:** 7FD-G011-08

**Oven Profile:**

**Carrier Gas:** Constant Flow Helium, 1 mL/min


**Injection:** Split 20.0000000000:1 0 µL @ 230°C

**Detection:** Refractive Index (330°C)

15867



Products used in this application:

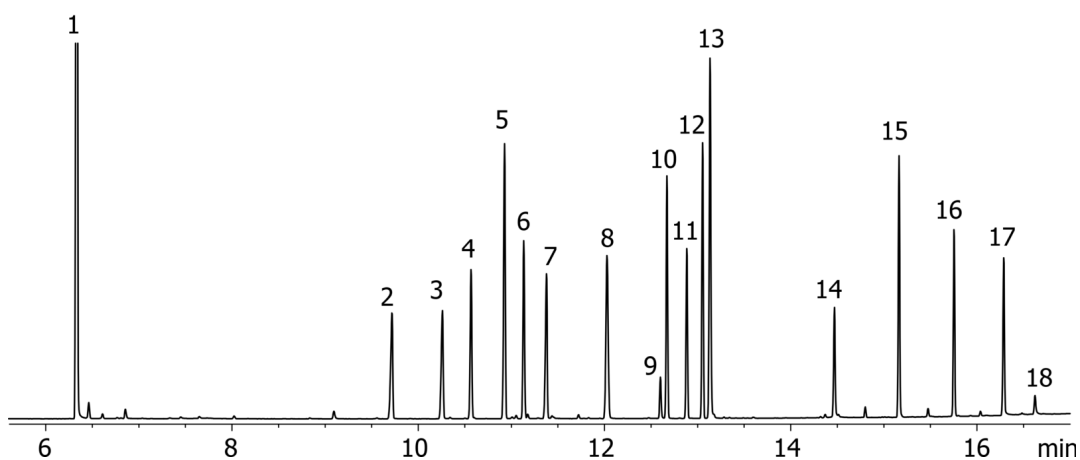

## Drug Screen on Zebron ZB-1MS 20m x 0.18mm x 0.18um

### ANALYTES:

- 1 Nicotine
- 2 Pentobarbital
- 3 Caffeine
- 4 Glutethimide
- 5 Lidocaine
- 6 Phencyclidine
- 7 Phenobarbital
- 8 EDDP
- 9 Methaqualone
- 10 Methadone
- 11 Cocaine
- 12 Trimipramine
- 13 Desipramine
- 14 Heroin
- 15 Fentanyl
- 16 Ibogaine
- 17 Triazolam
- 18 LSD

