

Drug Screen on Zebron ZB-1MS 12m x 0.20mm x 0.33um

Column: Zebron[™] ZB-1MS, GC Cap. Column 12 m x 0.20 mm x 0.33 µm, Ea

Phase: 100% Dimethylpolysiloxane

Dimensions: 12 meters x 0.2 mm x 0.33 µm

Order No: 7DE-G011-14

Oven Profile:

Carrier Gas: Constant Flow Helium, 1 mL/min


Injection: Split 20.000000000:1 0 µL @ 230°C

Detection: Refractive Index (330°C)

15869



Products used in this application:

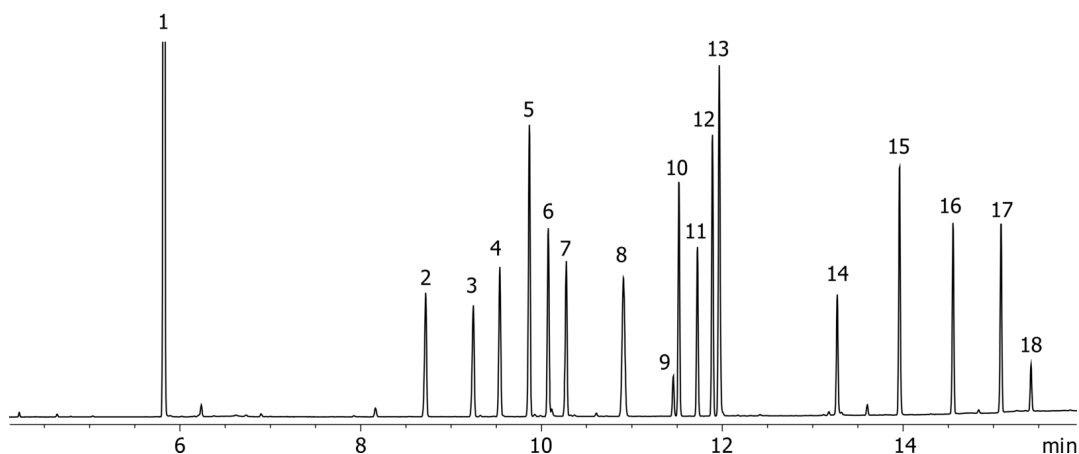

Drug Screen on Zebron ZB-1MS 12m x 0.20mm x 0.33um

ANALYTES:

- 1 Nicotine
- 2 Pentobarbital
- 3 Caffeine
- 4 Glutethimide
- 5 Lidocaine
- 6 Phencyclidine
- 7 Phenobarbital
- 8 EDDP
- 9 Methaqualone
- 10 Methadone
- 11 Cocaine
- 12 Trimipramine
- 13 Desipramine
- 14 Heroin
- 15 Fentanyl
- 16 Ibogaine
- 17 Triazolam
- 18 LSD

