


EPA Method 8081 (Organochlorine Pesticides) on ZB-1701P**Column:** Zebron™ ZB-1701P, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 14% Cyanopropylphenyl 86% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G012-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.6 mL/min**Injection:** Splitless 0.000000000:1 0 µL @ 240°C**Detection:** Refractive Index (300°C)**Analyst Note:** All analytes are 20ppm**Products used in this application:**

EPA Method 8081 (Organochlorine Pesticides) on ZB-1701P

ANALYTES:

- 1 1-Bromo-2-nitrobenzene
- 2 Tetrachloro-m-xylene
- 3 alpha-BHC
- 4 Pentachloronitrobenzene
- 5 gama-BHC
- 6 Heptachlor
- 7 Aldrin
- 8 b-BHC
- 9 delta-BHC
- 10 Heptachlor epoxide
- 11 Endosulfan I
- 12 gamma-Chlordane
- 13 alpha-Chlordane
- 14 4,4'-DDE
- 15 Dieldrin
- 16 Endrin
- 17 4,4'-DDD
- 18 Endosulfan II
- 19 4,4'-DDT
- 20 Endrin aldehyde
- 21 Methoxychlor
- 22 Endosulfan sulfate
- 23 Endrin ketone
- 24 Decachlorobiphenyl

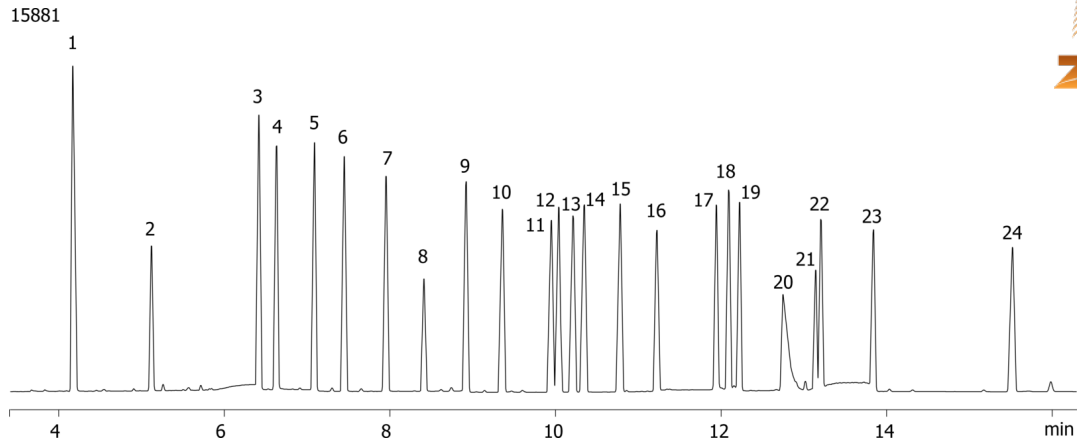