


**EPA Method 8081 - Organochlorine Pesticides on ZB-1ms****Column:** Zebron<sup>™</sup> ZB-1MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G011-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Splitless 0.000000000:1 0 µL @ 250°C**Detection:** Refractive Index (330°C)**Analyst Note:** All analytes are 20ppm.**Products used in this application:**

## EPA Method 8081 - Organochlorine Pesticides on ZB-1ms

### ANALYTES:

- 1 1-Bromo-2-Nitrobenzene
- 2 Tetrachloro-m-xylene
- 3 alpha-BHC
- 4 b-BHC
- 5 gama-BHC
- 6 delta-BHC
- 7 Pentachloronitrobenzene
- 8 Heptachlor
- 9 Aldrin
- 10 Heptachlor epoxide
- 11 gamma-Chlordane
- 12 alpha-Chlordane
- 13 Endosulfan I
- 14 4,4'-DDE
- 15 Dieldrin
- 16 Endrin
- 17 Endosulfan II
- 18 4,4'-DDD
- 19 Endrin aldehyde
- 20 Endosulfan sulfate
- 21 4,4'-DDT
- 22 Endrin ketone
- 23 Methoxychlor
- 24 Decachlorobiphenyl

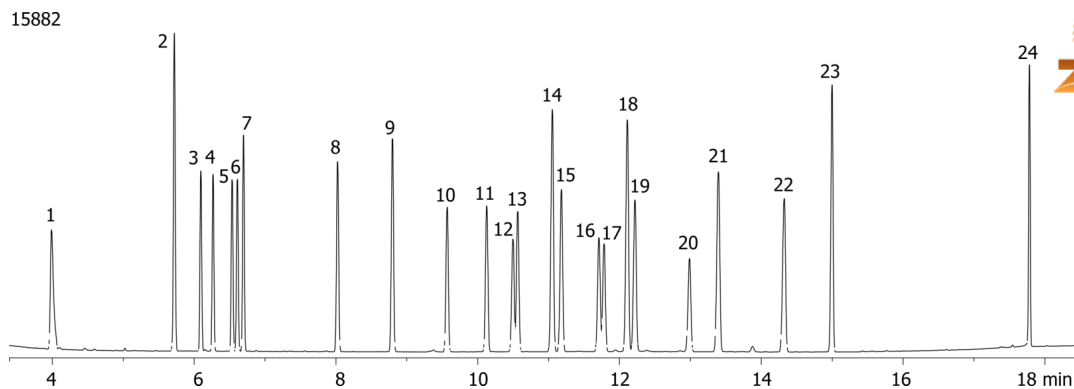