

Food Additives using Gemini 5u C6-Phenyl

Column: Gemini® 5 µm C6-Phenyl 110 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4444-E0**Elution Type:** Isocratic**Eluent A:** 20mM Potassium Phosphate pH 2.5**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	75	25

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Gemini C6-Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-7915

- Holder Part No.: KJ0-4282

15902



Products used in this application:

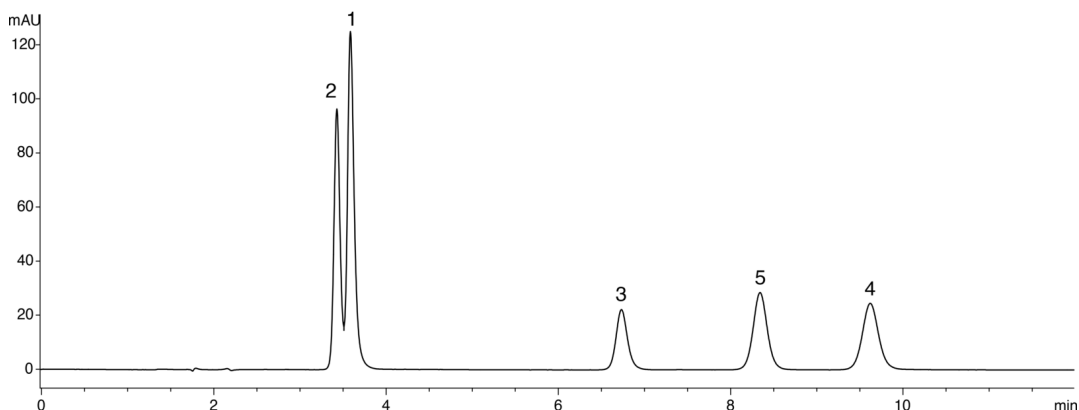

ANALYTES:

- 1 p-Hydroxybenzoic acid
- 2 Saccharine
- 3 Sorbic acid
- 4 Methyl paraben
- 5 Dehydroacetic acid

