


**Ethylene and Propylene Glycols on ZB-WAXplus****Column:** Zebtron<sup>™</sup> ZB-WAXplus<sup>™</sup>, GC Cap. Column 30 m x 0.25 mm x 0.25  $\mu$ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25  $\mu$ m**Order No:** 7HG-G013-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Splitless 0.0000000000:1 0  $\mu$ L @ 140°C**Detection:** Refractive Index (250°C)**Analyst Note:** Analytes are 100ppm in HPLC water**Products used in this application:****ANALYTES:**

- 1 Propylene glycol
- 2 Ethylene glycol
- 3 HPLC water contaminant

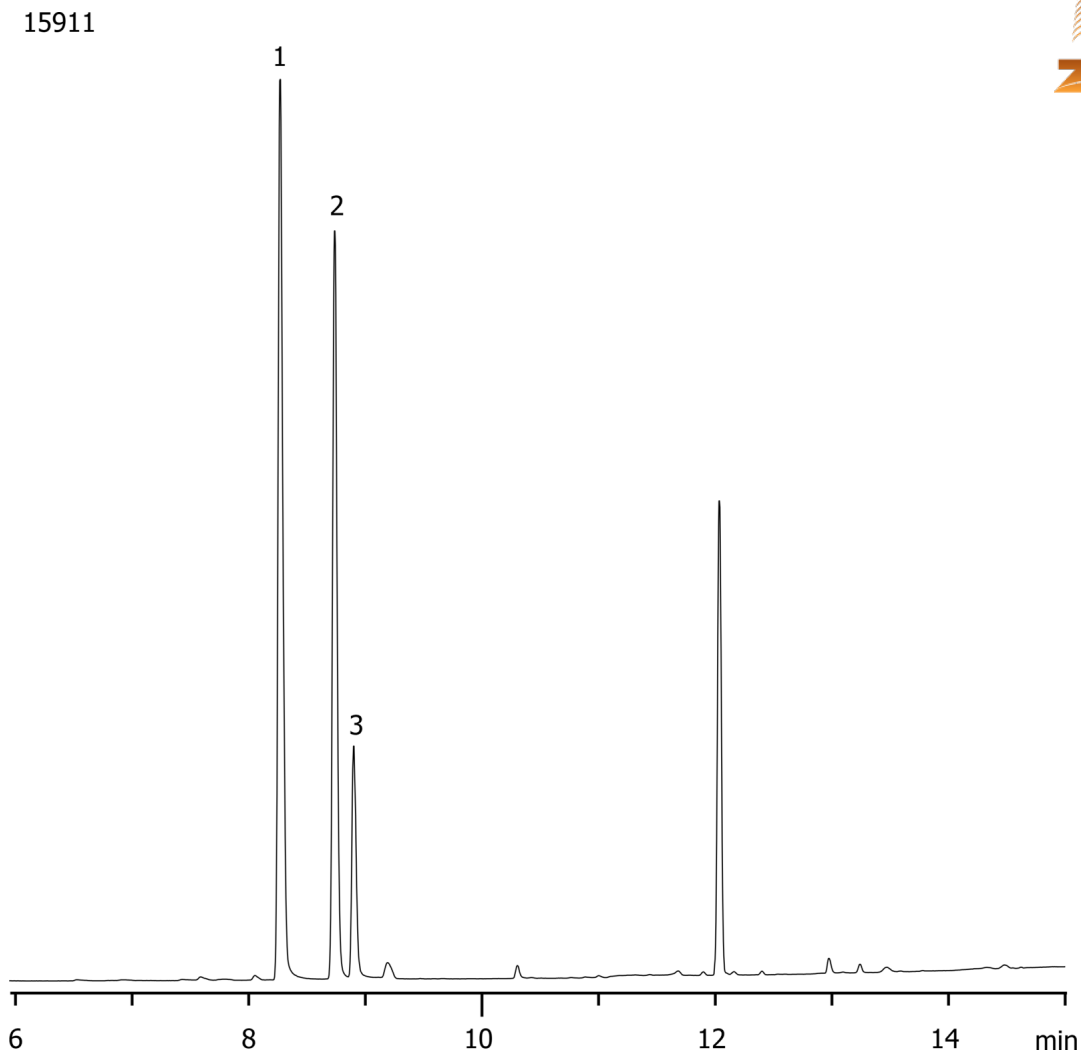