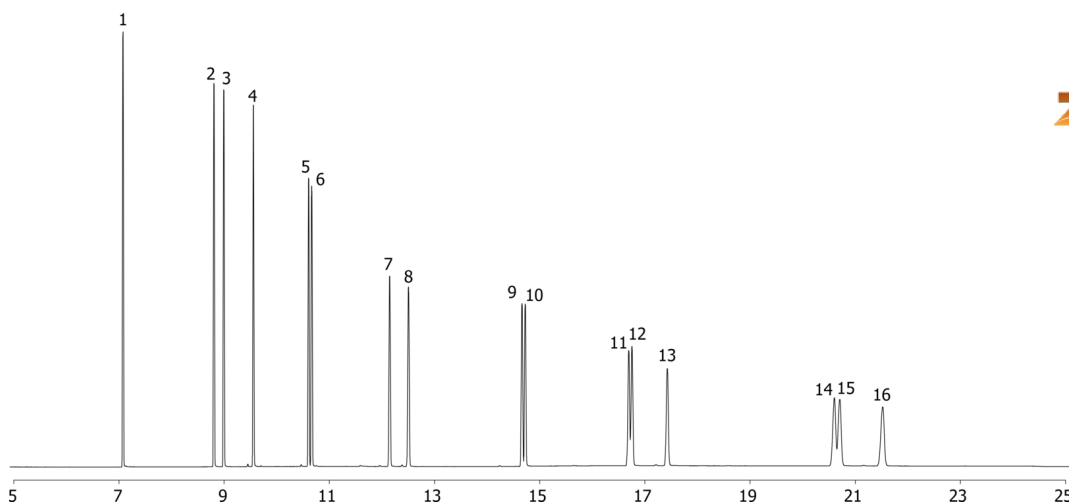


EPA 610 (Polynuclear Aromatic Hydrocarbons - PAHs) on ZB-5MS**Column:** Zebtron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:****Carrier Gas:** Constant Pressure Hydrogen, 1.2 mL/min**Injection:** Split 125.0000000000:1 0 µL @ 140°C**Detection:** Refractive Index (325°C)

15914

**Products used in this application:**

EPA 610 (Polynuclear Aromatic Hydrocarbons - PAHs) on ZB-5MS

ANALYTES:

- 1 Naphthalene
- 2 Acenaphthylene
- 3 Acenaphthene
- 4 Fluorene
- 5 Phenanthrene
- 6 Anthracene
- 7 Fluoranthene
- 8 Pyrene
- 9 Benzo[a]anthracene
- 10 Chrysene
- 11 Benzo[b]fluoranthene
- 12 Benzo[k]fluoranthene
- 13 Benzo[a]pyrene
- 14 Indeno[1,2,3-cd]pyrene
- 15 Dibenzo[a,h]anthracene
- 16 Benzo[g,h,i]perylene

