

## Basic Compounds at Low pH Using Luna 5u C18(2) Axia Packed 50x21.2mm

**Column:** Luna® 5 µm C18(2) 100 Å, LC Column 50 x 21.2 mm, AXIA™ Packed, Ea

**Dimensions:** 50 x 21.2 mm ID

**Order No:** 00B-4252-P0-AX

**Elution Type:** Gradient

**Eluent A:** 0.5% TFA in Water

**Eluent B:** 0.5% TFA in Acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 95    | 5     |
|                   | 2        | 7          | 5     | 95    |

**Flow Rate:** 30 mL/min

**Col. Temp.:** ambient

**Detection:** Refractive Index @ 0.0000000000 (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard PREP Cartridge, C18 15 x 21.2mm ID, Ea Part No.: AJ0-7839


- Holder Part No.: AJ0-8223



Products used in this application:



15929

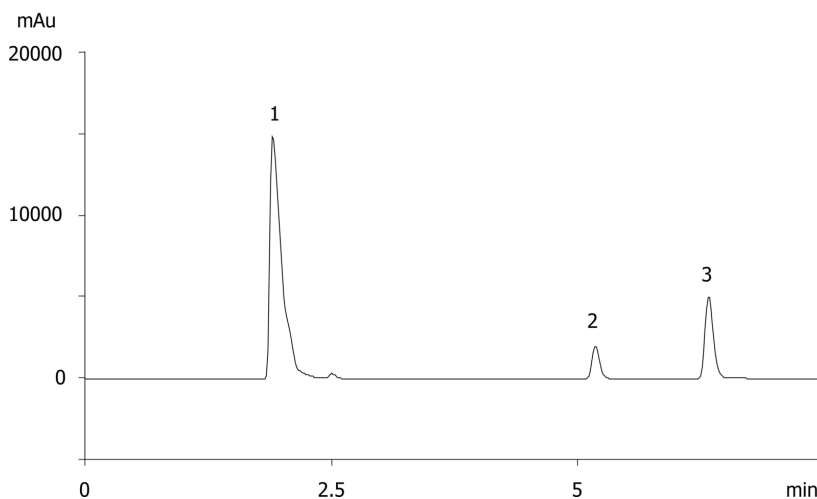

### ANALYTES:

- 1 Triprolidine
- 2 Methacycline
- 3 Amitriptyline

