


Common Solvents on ZB-WAX**Column:** ZebronTM ZB-WAX, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G007-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 25.0000000000:1 0 μ L @ 100°C**Detection:** Refractive Index (0°C)**Analyst Note:** 4 μ L Headspace split 25:1 at 100°C**Products used in this application:**

Common Solvents on ZB-WAX

ANALYTES:

- 1 Acetone
- 2 Ethyl acetate
- 3 Isopropyl acetate
- 4 Methanol
- 5 Methyl ethyl ketone
- 6 Isopropyl alcohol
- 7 Ethanol
- 8 n-Propyl acetate
- 9 Methyl isobutyl ketone
- 10 n-Propanol
- 11 Toluene
- 12 n-Butyl acetate
- 13 Methoxypropanol (Propylene glycol monomethyl ether)
- 14 3-Ethoxy propanol (Propylene glycol Ä-monoethyl ether)

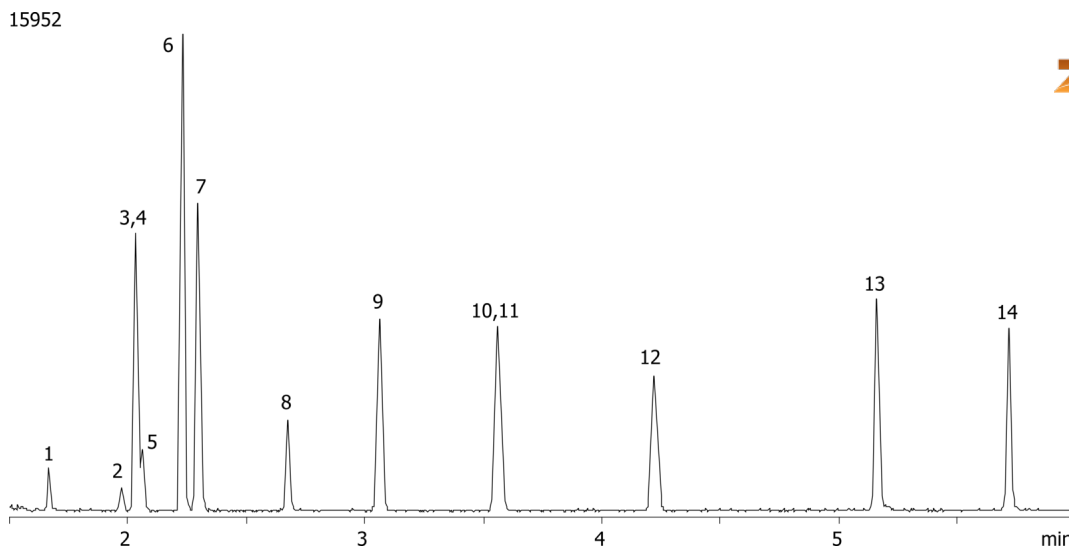