


Common Drug Screen on ZB-5MS**Column:** Zebtron[™] ZB-5MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 20 meters x 0.18 mm x 0.18 µm**Order No:** 7FD-G010-08**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 225°C**Detection:** Refractive Index (250°C)**Analyst Note:** All analytes are 25ppm except nicotine at 100ppm**Products used in this application:**

Common Drug Screen on ZB-5MS

ANALYTES:

- 1 Nicotine
- 2 Acetophenone
- 3 Ibuprofen
- 4 Allobarbitol
- 5 Aprobarbitol
- 6 Butalbitol
- 7 Phenacetin
- 8 Amobarbitol
- 9 Pentobarbitol
- 10 Meprobamate
- 11 Secobarbitol
- 12 Caffeine
- 13 Benzphetamine
- 14 Hexobarbitol
- 15 Dimenhydrinate
- 16 Doxylamine
- 17 Phenobarbitol
- 18 Methapyrilene
- 19 Chlorpheniramine
- 20 Procaine
- 21 Brompheniramine
- 22 Cocaine
- 23 Chlorcyclizine
- 24 Benzactyzine
- 25 Codeine
- 26 Diazepam
- 27 Morphine
- 28 Hydrocodone
- 29 Oxymorphone
- 30 Heroin

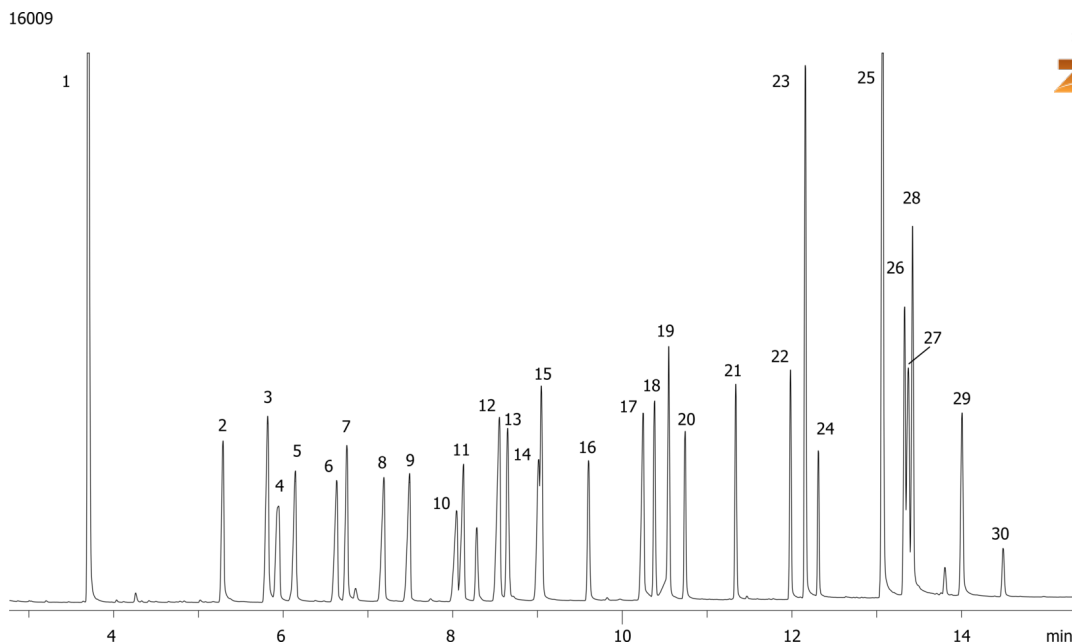