


**Common Drug Screen on ZB-50****Column:** Zebron™ ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 50% Phenyl 50% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G004-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2.2 mL/min**Injection:** Split 10.000000000:1 0 µL @ 275°C**Detection:** Refractive Index (300°C)**Analyst Note:** All analytes are 25ppm except nicotine at 100ppm**Products used in this application:**

## Common Drug Screen on ZB-50

### ANALYTES:

- 1 Nicotine
- 2 Ibuprofen
- 3 Allobarbitol
- 4 Acetomenophen
- 5 Aprobital
- 6 Butalital
- 7 Amobarbitol
- 8 Pentobarbitol
- 9 Phenacetin
- 10 Benzphetamine
- 11 Secobarbitol
- 12 Dimenhydrinate
- 13 Meprobamate
- 14 Doxylamine
- 15 Hexobarbitol
- 16 Caffeine
- 17 Chlorpheniramine
- 18 Methapyrilene
- 19 Phenobarbitol
- 20 Procaine
- 21 Brompheniramine
- 22 Chlorcyclizine
- 23 Cocaine
- 24 Benactyzine
- 25 Codeine
- 26 Diazepam
- 27 Morphine
- 28 Hydrocodone
- 29 Oxymorphone
- 30 Heroin

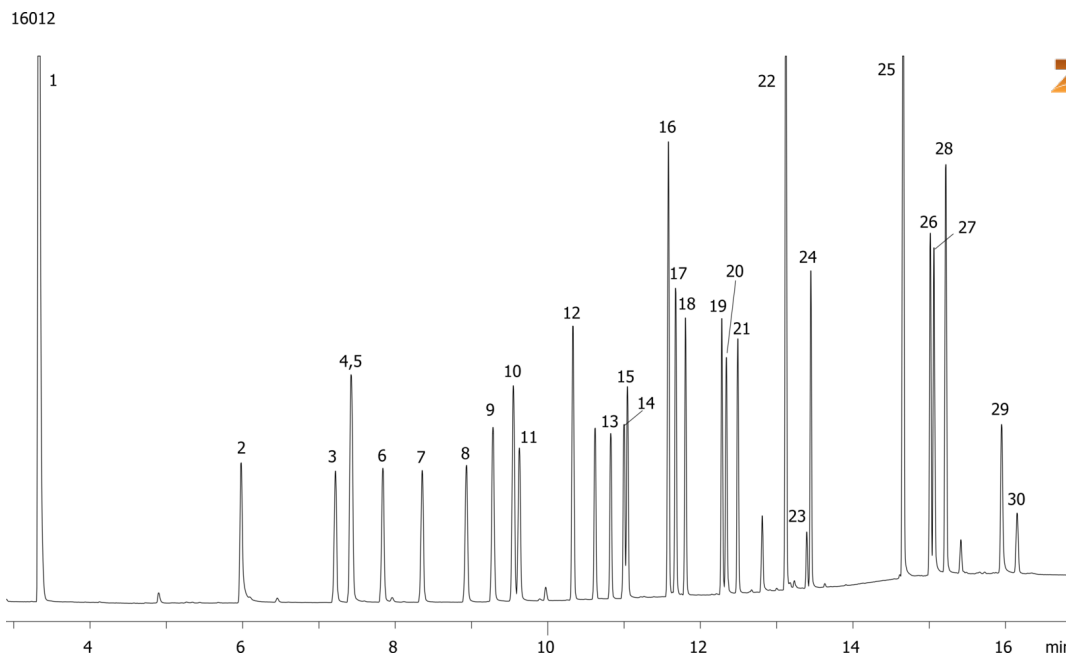