


**Peppermint Oil on ZB-1MS****Column:** Zebtron<sup>™</sup> ZB-1MS, GC Cap. Column 10 m x 0.10 mm x 0.10 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.1 mm x 0.1 µm**Order No:** 7CB-G011-02**Oven Profile:****Carrier Gas:** Constant Flow Helium, 0.3 mL/min**Injection:** Split 120.0000000000:1 0 µL @ 160°C**Detection:** Refractive Index (0°C)**Analyst Note:** 10% Peppermint Oil in Dichloromethane.**Products used in this application:**

16018

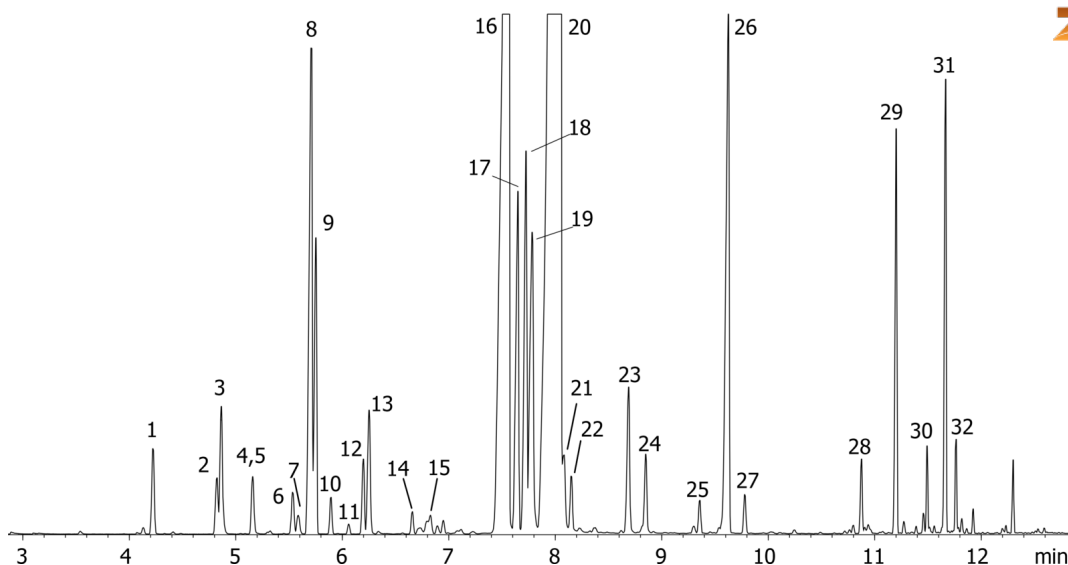

## Peppermint Oil on ZB-1MS

### ANALYTES:

- 1 alpha-Pinene
- 2 Sabinene
- 3 beta-Pinene
- 4 beta-Myrcene
- 5 3-Octanol
- 6 a-Terpenene
- 7 Cymene
- 8 Eucalyptol
- 9 D-Limonene
- 10 cis-Ocimene
- 11 trans-Ocimene
- 12 g-Terpinene
- 13 cis-Sabinene hydrate
- 14 b-Terpineol
- 15 Linalool
- 16 Menthone
- 17 Isomenthone
- 18 Menthonefuran
- 19 Neomenthol
- 20 Menthol
- 21 Neoisomenthol
- 22 a-Terpineol
- 23 Pulegone
- 24 Piperitone
- 25 Neomenthyl acetate
- 26 Menthyl acetate
- 27 Isomenthyl acetate
- 28 Bourbonene
- 29 Caryophyllene
- 30 Farnesene
- 31 Germacrene
- 32 Elemene

