


European Red List Volatiles on ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2.7 mL/min**Injection:** Split 50.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (260°C)**Analyst Note:** Static Headspace: 8µL of 12.5% volatiles mix sample in a 0.18mL limited volume insert inside a 2.0mL vial.**Products used in this application:****ANALYTES:**

- 1** 1,1-Dichloroethylene
- 2** 1,1-Dichloroethane
- 3** Chloroform
- 4** 1,1,1-Trichloroethane
- 5** 1,2-Dichloroethane
- 6** Carbon tetrachloride
- 7** Trichloroethylene
- 8** Tetrachloroethylene

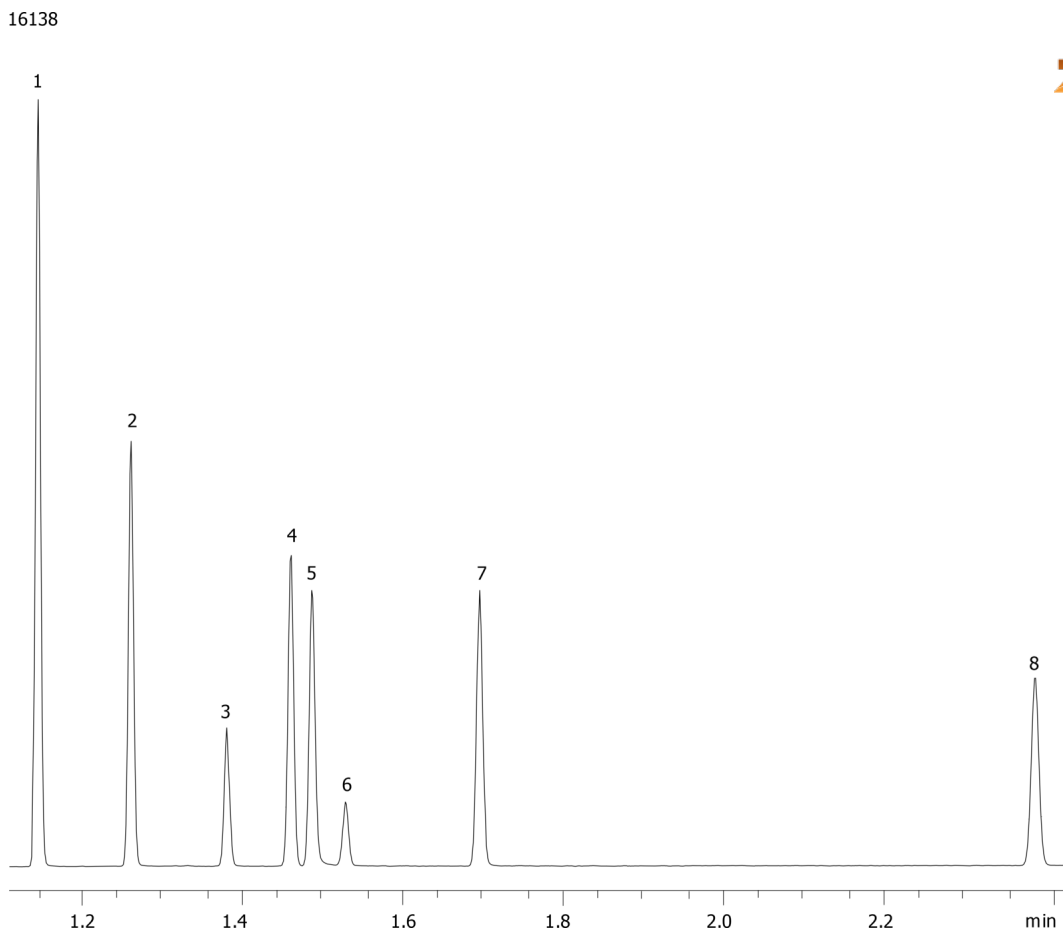