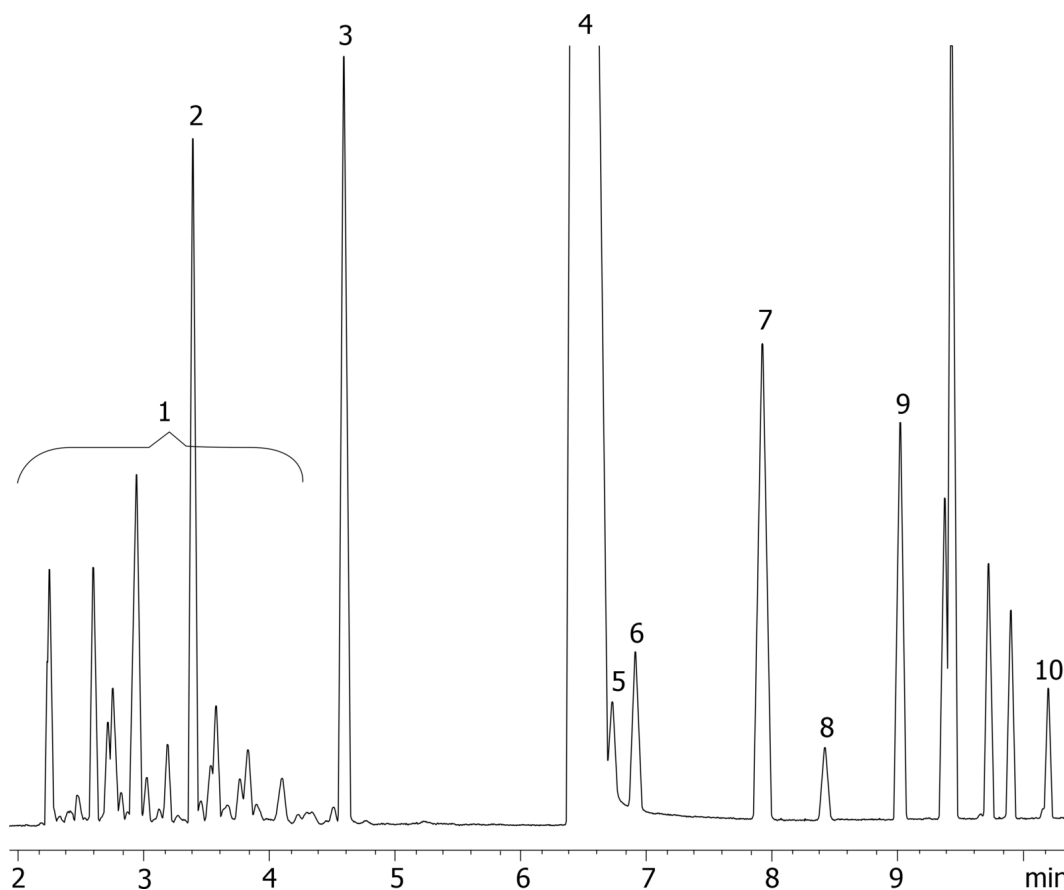


ASTM Method D5060: Determining Impurities in High Purity Ethylbenzene on ZB-FFAP**Column:** Zebron™ ZB-FFAP, GC Cap. Column 50 m x 0.32 mm x 0.50 µm, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 50 meters x 0.32 mm x 0.5 µm**Order No:** 7JM-G009-17**Oven Profile:****Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** Split 200.0000000000:1 0 µL @ 160°C**Detection:** Refractive Index (240°C)

16145

**Products used in this application:**

ASTM Method D5060: Determining Impurities in High Purity Ethylbenzene on ZB-FFAP

ANALYTES:

- 1 Aliphatic hydrocarbons
- 2 Benzene
- 3 Toluene
- 4 Ethylbenzene
- 5 p-Xylene
- 6 m-Xylene
- 7 Cumene
- 8 o-Xylene
- 9 Propylbenzene
- 10 Styrene

