


DDD, DDT, And DDE Isomers on ZB-MR1**Column:** Zebtron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G016-11**Oven Profile:****Carrier Gas:** Constant Pressure Helium, 0.8 mL/min**Injection:** Splitless 0.0000000000:1 0 µL @ 260°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes are 25 ng on column**Products used in this application:****ANALYTES:**

- 1 o,p'-DDE
- 2 p,p'-DDE
- 3 o,p'-DDD
- 4 o,p'-DDT
- 5 p,p'-DDD
- 6 p,p'-DDT

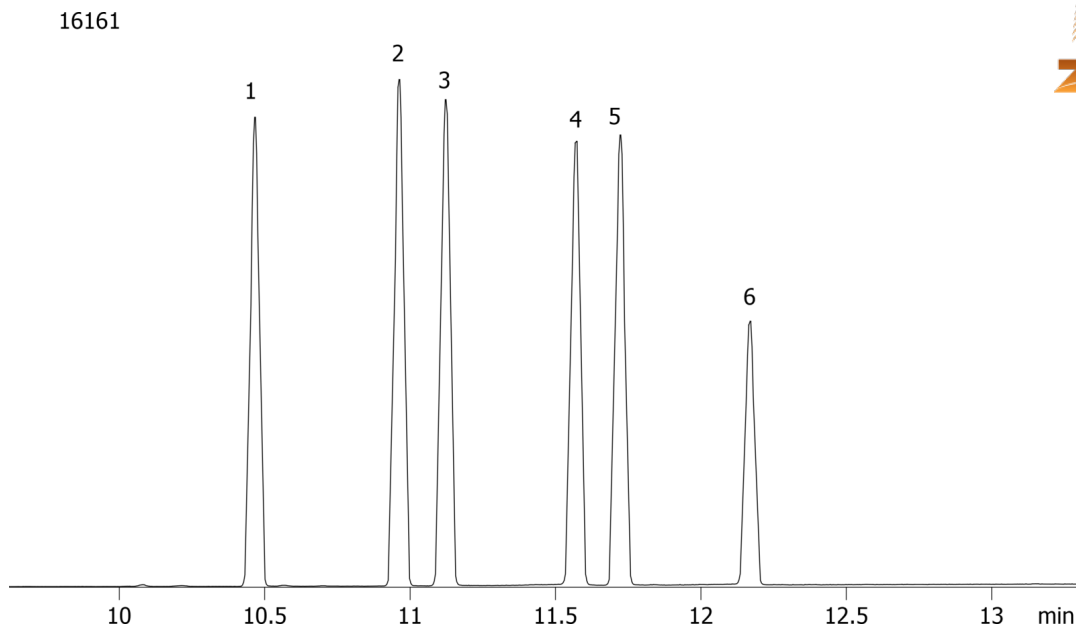