


DDD, DDT, And DDE Isomers on ZB-MR2**Column:** Zebron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:****Carrier Gas:** Constant Pressure Helium, 0.8 mL/min**Injection:** Splitless 0.000000000:1.0 μ L @ 260°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes are 25 ng on column**Products used in this application:****ANALYTES:**

- 1** o,p'-DDE
- 2** p,p'-DDE
- 3** o,p'-DDD
- 4** o,p'-DDT
- 5** p,p'-DDD
- 6** p,p'-DDT

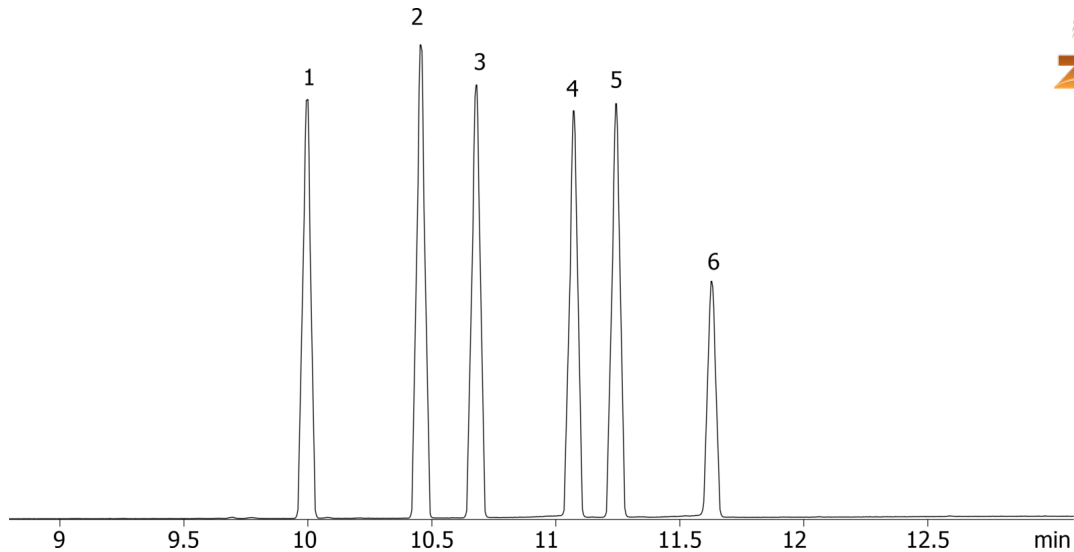