


EPA 8141: Organophosphorous Pesticides on ZB-MR2**Column:** Zebron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G017-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2.8 mL/min**Injection:** On-Column 0.000000000:1 0 µL @ 103°C**Detection:** Refractive Index (340°C)**Analyst Note:** Analytes are 2ppm in Dichloromethane. Columns connected via a 5 m guard and 'Y' connector.**Products used in this application:**

EPA 8141: Organophosphorous Pesticides on ZB-MR2**ANALYTES:**

- | | |
|-------------------------------|---------------------------|
| 1 Dichlorvos | 41 Carbophenothion |
| 2 Mevinphos | 42 Famphur |
| 3 Trichlorfon | 43 EPN |
| 4 TEPP | 44 Phosmet |
| 5 Demeton | 45 Leptophos |
| 6 Thionazin | 46 Azinphos methyl |
| 7 Ethoprop | 47 Azinphos ethyl |
| 8 Sulfotep | 48 Coumaphos |
| 9 Naled | |
| 10 Phorate | |
| 11 Dicrotophos | |
| 12 Demeton | |
| 13 Monocrotophos | |
| 14 Terbufos | |
| 15 Diazinon | |
| 16 Fonofos | |
| 17 Dimethoate | |
| 18 Disulfoton | |
| 19 Phosphamidon | |
| 20 Dichlofenthion | |
| 21 Phosphamidon | |
| 22 Chlorpyrifos methyl | |
| 23 Ronnel | |
| 24 Aspon | |
| 25 Methyl parathion | |
| 26 Malathion | |
| 27 Fenitrothion | |
| 28 Chlorpyrifos | |
| 29 Trichloronate | |
| 30 Fenthion | |
| 31 Parathion | |
| 32 Merphos | |
| 33 Chlorfenvinphos | |
| 34 Crotoxyphos | |
| 35 Tetrachlorvinphos | |
| 36 Tokuthion | |
| 37 Merphos oxide | |
| 38 Ethion | |
| 39 Fensulfothion | |
| 40 contaminant | |

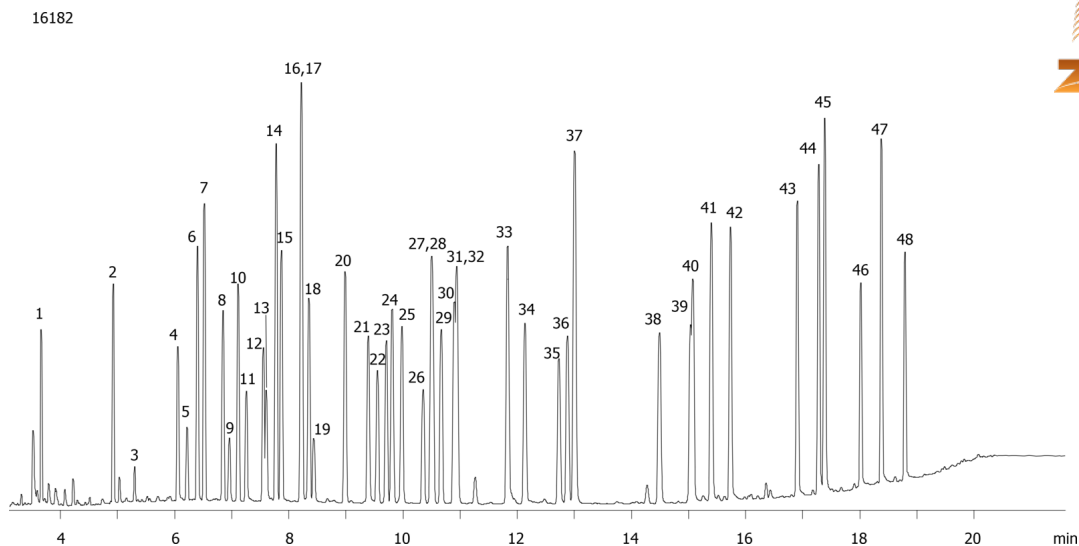