


Drug Screen on ZB-1HT**Column:** Zebron[™] ZB-1HT, GC Cap. Column 30 m x 0.25 mm x 0.10 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.1 µm**Order No:** 7HG-G014-02**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.1 mL/min**Injection:** Split 15.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (310°C)**Analyst Note:** Analytes were 25ppm.**Products used in this application:**

Drug Screen on ZB-1HT

ANALYTES:

- 1 Nicotine
- 2 contaminant
- 3 Acetaminophen
- 4 Allobarbitol
- 5 Aprobarbitol
- 6 Ibuprofen
- 7 Butabarbital
- 8 Phenacetin
- 9 Amobarbital
- 10 Pentobarbital
- 11 Meprobamate
- 12 Secobarbital
- 13 Caffeine
- 14 Methyl benzilate
- 15 Benzphetamine
- 16 Hexobarbital
- 17 Dimenhydrinate
- 18 Doxylamine
- 19 Phenobarbital
- 20 Methapyrilene
- 21 Chlorpheniramine
- 22 Brompheniramine
- 23 Cocaine
- 24 Chlorcyclizine
- 25 Codeine
- 26 Diazepam
- 27 Hydrocodone
- 28 Oxymorphone

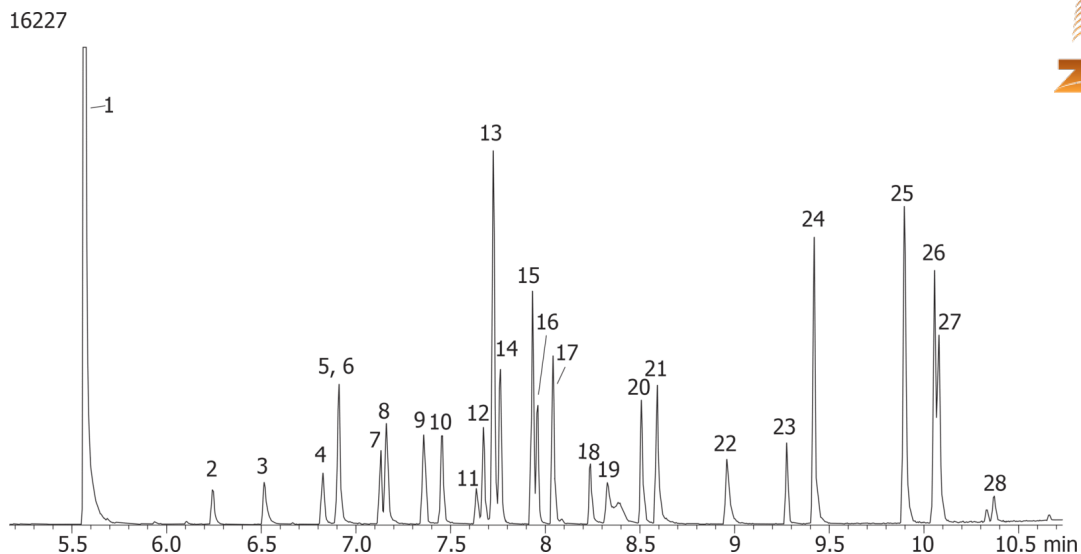