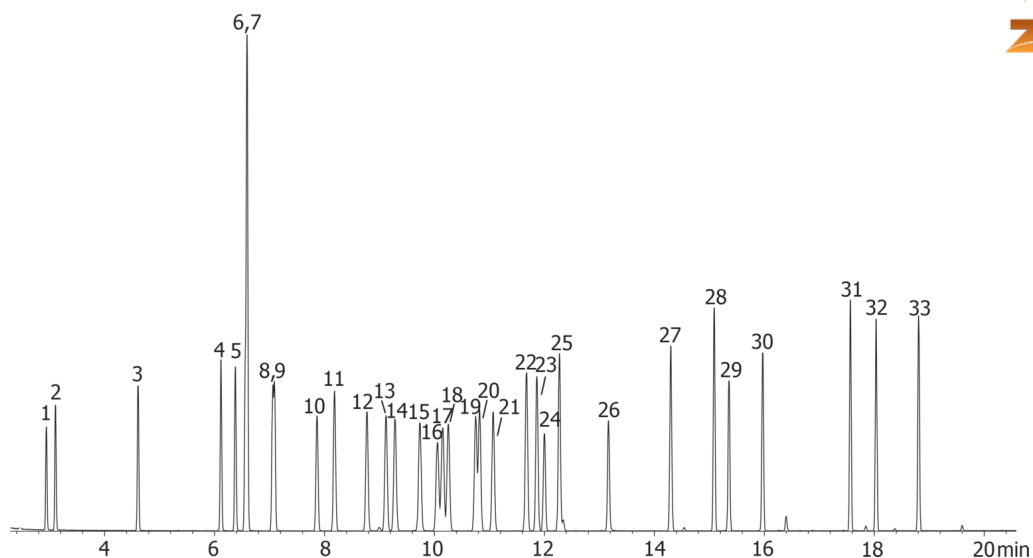


Aromatics on ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.32 mm x 1.00 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 1 µm**Order No:** 7HM-G010-22**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.6 mL/min**Injection:** Split 20.000000000:1 0 µL @ 250°C**Detection:** Refractive Index (275°C)**Analyst Note:** The analytes are 100ppm in Methylene Chloride

16229

**Products used in this application:**

Aromatics on ZB-5MS

ANALYTES:

- 1 Benzene
- 2 Fluorobenzene
- 3 Toluene
- 4 Chlorobenzene
- 5 Ethylbenzene
- 6 m-Xylene
- 7 p-Xylene
- 8 Styrene
- 9 o-Xylene
- 10 1-Methylethylbenzene
- 11 Bromobenzene
- 12 Propylbenzene
- 13 4-Ethyltoluene
- 14 1,3,5-Trimethylbenzene
- 15 a-Methylstyrene
- 16 Ethoxybenzene
- 17 tert-Butylbenzene
- 18 1,2,4-Trimethylbenzene
- 19 2-Methylpropylbenzene
- 20 1-Methylpropylbenzene
- 21 1,4-Dichlorobenzene
- 22 1,2-Dichlorobenzene
- 23 Iodobenzene
- 24 Diethylbenzene
- 25 Butylbenzene
- 26 Nitrobenzene
- 27 1,3,5-Trichlorobenzene
- 28 1,4-Diisopropylbenzene
- 29 2-Nitro-m-xylene
- 30 1,2,3-Trichlorobenzene
- 31 1-Bromo-3-phenyl-propane
- 32 1,2,4,5-Tetrachlorobenzene
- 33 1,2,3,4-Tetrachlorobenzene

