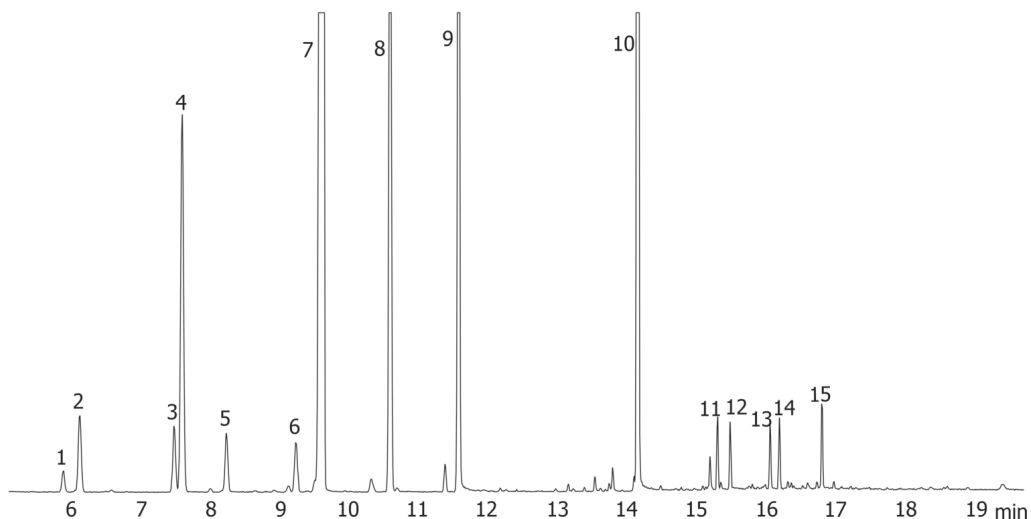


Bergamot Oil on ZB-1**Column:** Zebtron™ ZB-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G001-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** On-Column 0.000000000:1 0 µL @ 63°C**Detection:** Refractive Index (240°C)**Analyst Note:** The solution analyzed was 1000ppm Bergamot oil in Ethyl acetate

16230

**Products used in this application:**

Bergamot Oill on ZB-1

ANALYTES:

- 1 a-Thujene
- 2 a-Pinene
- 3 Sabinene
- 4 beta-Pinene
- 5 beta-Myrcene
- 6 p-Cymene
- 7 Limonene
- 8 g-Terpinene
- 9 Linalool acetate
- 10 Linalyl acetate
- 11 Neryl acetate
- 12 Geranyl acetate
- 13 b-Caryophyllene
- 14 alpha-Bergamotene
- 15 b-Bisabolene

