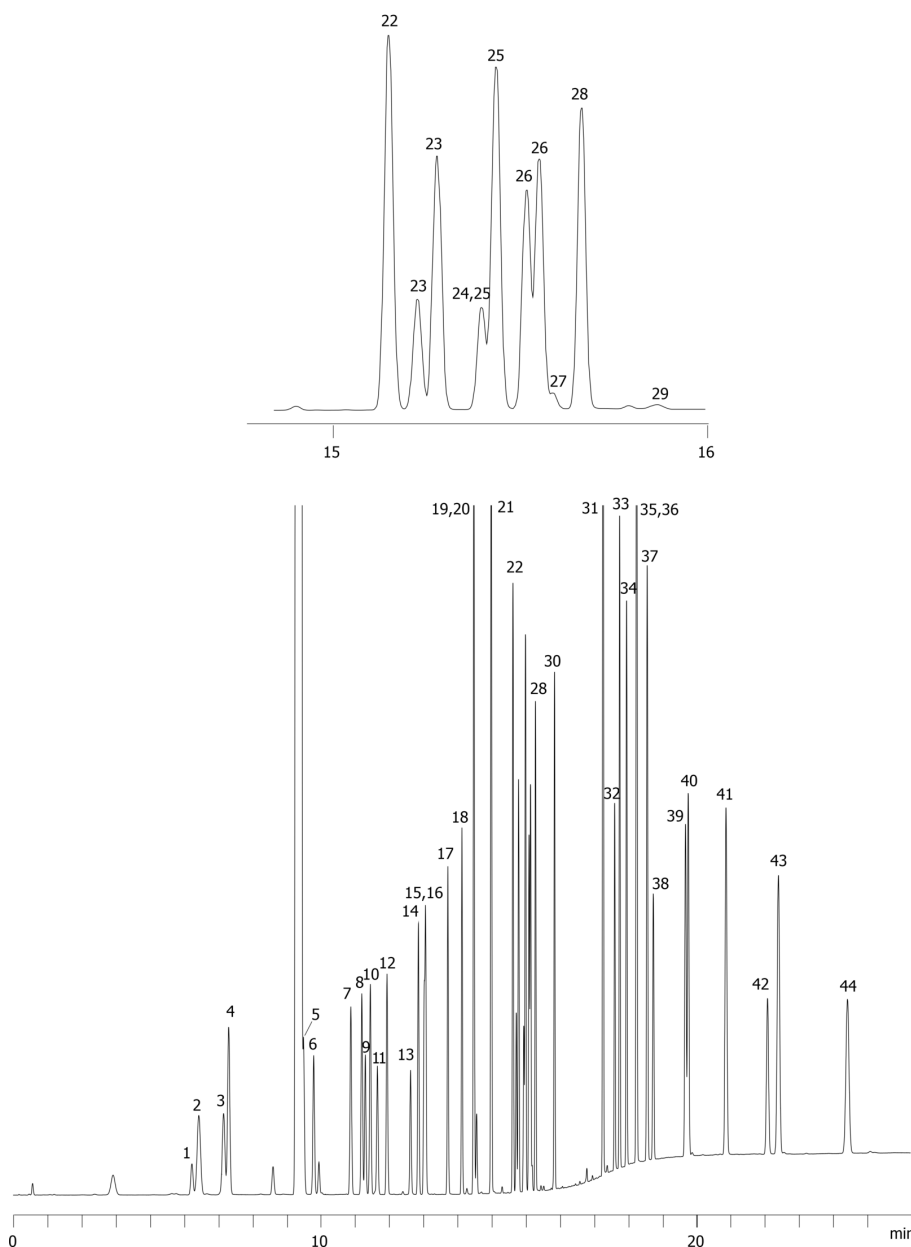


Industrial Chemicals: Alcohols on ZB-WAXplus**Column:** Zebron™ ZB-WAXplus™, GC Cap. Column 30 m x 0.53 mm x 1.00 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.53 mm x 1 µm**Order No:** 7HK-G013-22**Oven Profile:****Carrier Gas:** Constant Flow Helium, 4.1 mL/min**Injection:** Split 10.000000000:1 0 µL @ 200°C**Detection:** Refractive Index (220°C)**Analyst Note:** Analytes are 100ppm in Carbon disulfide

16234

**Products used in this application:**

Industrial Chemicals: Alcohols on ZB-WAXplus**ANALYTES:**

1	Methanol	41	1-Tetradecanol
2	tert-Butanol	42	1,7-Heptanediol
3	Isopropanol	43	1-Pentadecanol
4	Ethanol	44	1-Hexadecanol
5	2-Butanol		
6	1-Propanol		
7	2-Methylpropanol		
8	1-Pentanol		
9	Allyl alcohol		
10	2-Pentanol		
11	1-Methoxy-2-propanol		
12	1-Butanol		
13	2-Methoxyethanol		
14	3-Hexanol		
15	Methallyl alcohol		
16	3-Methyl butanol		
17	Pentanol		
18	4-Heptanol		
19	2-Ethyl-1-butanol		
20	Cyclopentanol		
21	Hexanol		
22	Cyclohexanol		
23	2-Methylcyclohexanol		
24	1-Heptanol		
25	3-Methylcyclohexanol		
26	4-Methylcyclohexanol		
27	Glycidol		
28	2-Ethyl-1-hexanol		
29	Mercaptoethanol		
30	Octanol		
31	Decanol		
32	2-(2-Butoxyethoxy)ethanol		
33	1-Phenylethanol		
34	1-Undecanol		
35	Benzyl alcohol		
36	Phenylmethanol		
37	2-Phenylethanol		
38	1-Dodecanol		
39	1-Tridecanol		
40	3-Phenyl-1-propanol		

