


EPA 8081 including PCB congeners and o,p'-DDD on ZB-MR1**Column:** ZebtronTM ZB-MultiResidueTM-1, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G016-17**Oven Profile:****Carrier Gas:** Constant Pressure Helium, 3.6 mL/min**Injection:** On-Column 0.000000000:1 0 μ L @ 103°C**Detection:** Refractive Index (340°C)**Analyst Note:** Analytes are 50ppb in iso-octane**Products used in this application:**

EPA 8081 including PCB congeners and o,p'-DDD on ZB-MR1

ANALYTES:

- 1 alpha-BHC
- 2 gama-BHC
- 3 b-BHC
- 4 delta-BHC
- 5 Heptachlor
- 6 BZ# 52
- 7 Aldrin
- 8 Heptachlor epoxide
- 9 gamma-Chlordane
- 10 alpha-Chlordane
- 11 Endosulfan I
- 12 4,4'-DDE
- 13 2,4'-DDD
- 14 Dieldrin
- 15 BZ# 118
- 16 Endrin
- 17 4,4'-DDD
- 18 Endosulfan II
- 19 Endrin aldehyde
- 20 4,4'-DDT
- 21 Endosulfan sulfate
- 22 Methoxychlor
- 23 Endrin ketone

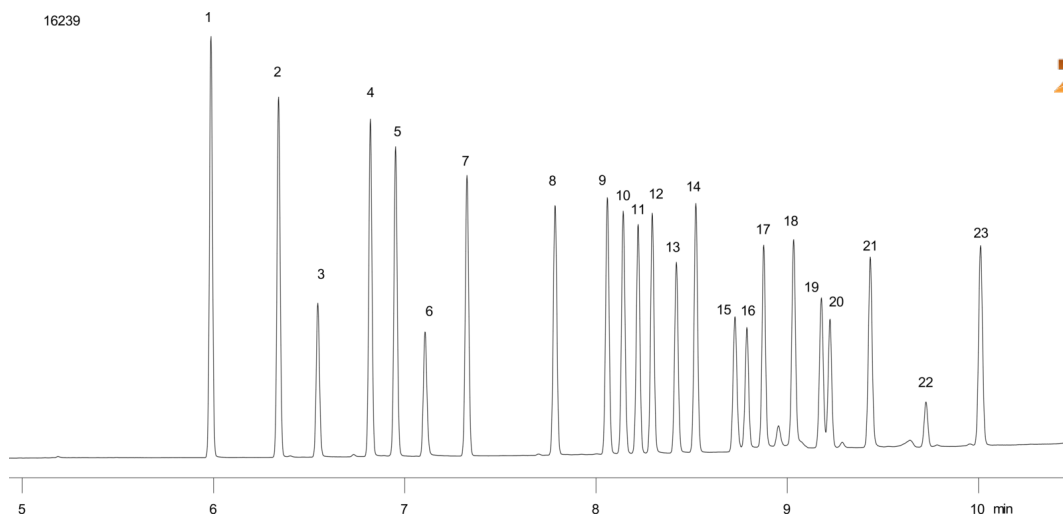