


ASTM Method D2887-97: Boiling Range Distribution of Petroleum Fractions on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** Split 10.000000000:1 0 μ L @ 300°C**Detection:** Refractive Index (350°C)**Analyst Note:** Diesel Fuel at 600 ppm in Hexane**Products used in this application:**

16255

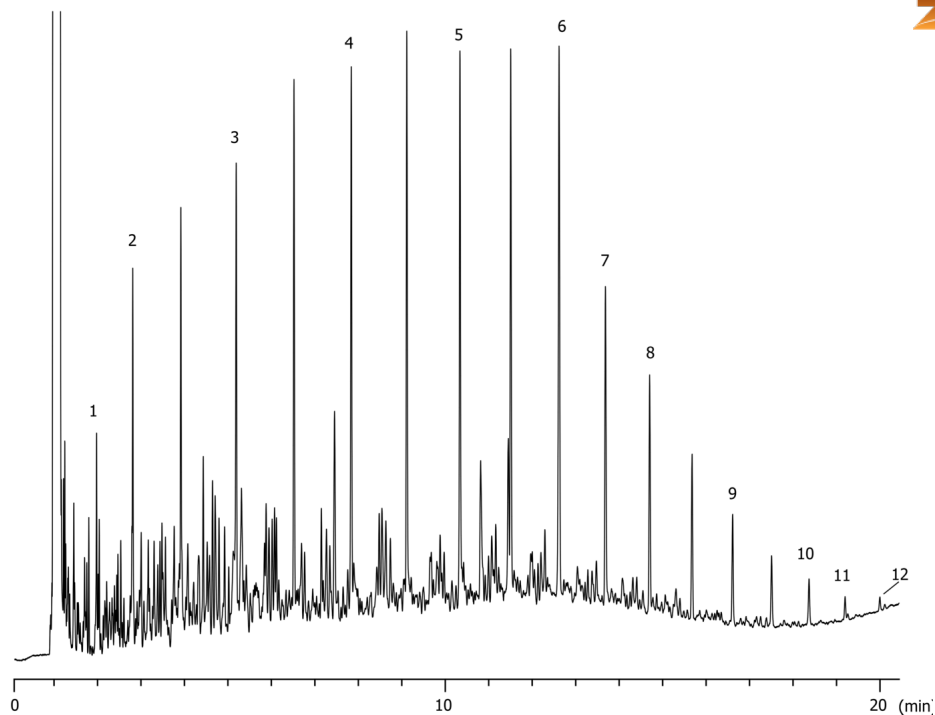

ASTM Method D2887-97: Boiling Range Distribution of Petroleum Fractions on ZB-MR1

ANALYTES:

- 1 n-Nonane (C9)
- 2 n-Decane (C10)
- 3 n-Dodecane (C12)
- 4 n-Tetradecane (C14)
- 5 n-Hexadecane (C16)
- 6 n-Octadecane (C18)
- 7 n-Nonadecane (C19)
- 8 n-Eicosane (C20)
- 9 n-Docosane (C22)
- 10 n-Tetracosane (C24)
- 11 n-Hexacosane (C26)
- 12 n-Octacosane (C28)

