


Vitamin E and Sterols on ZB-5**Column:** Zebtron[™] ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.10 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.1 µm**Order No:** 7HG-G002-02**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.8 mL/min**Injection:** Splitless 0.000000000:1 0 µL @ 220°C**Detection:** Refractive Index (340°C)**Analyst Note:** Recommended Accessories:
Recommended Liner: Zebtron PLUS Single Taper Z-Liner, 4 mm ID
Liner Part No.: AG2-0A13-05 (for Agilent systems)**Products used in this application:**

16331

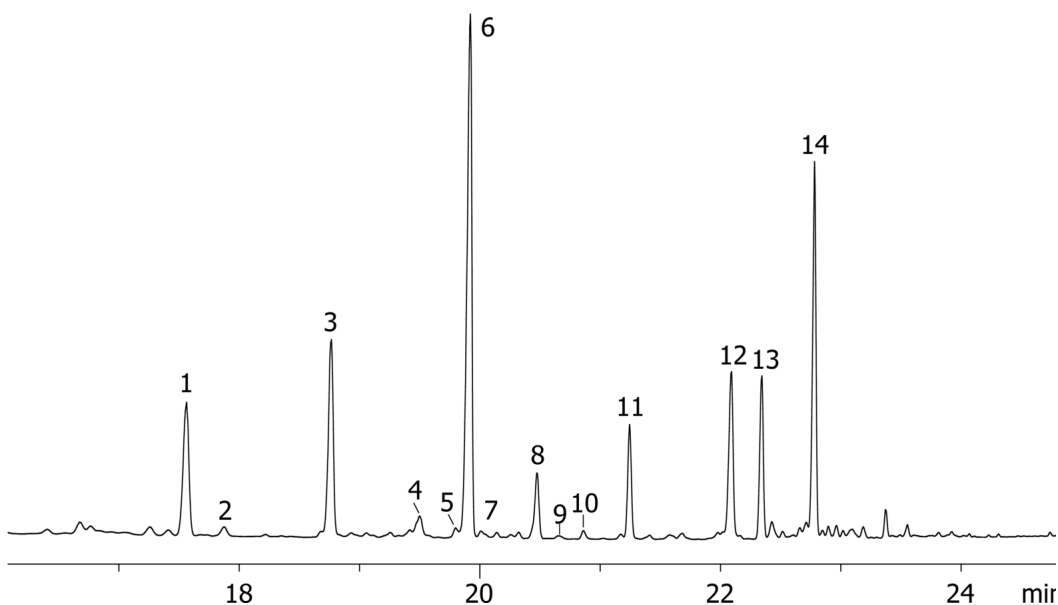

Vitamin E and Sterols on ZB-5

ANALYTES:

- 1 Squalene
- 2 Lignoceric acid
- 3 delta-Tocopherol
- 4 d-Tocomonoenol
- 5 Campesta-3,5-diene
- 6 gamma-Tocopherol
- 7 Stigmasta-3,5,22-triene
- 8 g-Tocomonoenol
- 9 Stigmasta-3,5-diene
- 10 Cholesterol
- 11 alpha-Tocopherol
- 12 Campesterol
- 13 Stigmasterol
- 14 b-Sitosterol

