


USP 467 Procedure A - Class 2 Mix B on ZB-624**Column:** Zebron[™] ZB-624, GC Cap. Column 30 m x 0.32 mm x 1.80 μ m, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 1.8 μ m**Order No:** 7HM-G005-31**Oven Profile:****Carrier Gas:** Constant Flow Helium, 35 cm/sec**Injection:** Headspace 5.000000000:1 0 μ L @ 140°C**Detection:** Refractive Index (250°C)**Analyst Note:** Sample Prepared as per USP Guidelines**Products used in this application:****ANALYTES:**

- 1 Hexane
- 2 Nitromethane
- 3 Chloroform
- 4 1,2-Dimethoxyethane
- 5 Trichloroethylene
- 6 Pyridine
- 7 Methyl butyl ketone
- 8 Tetralin

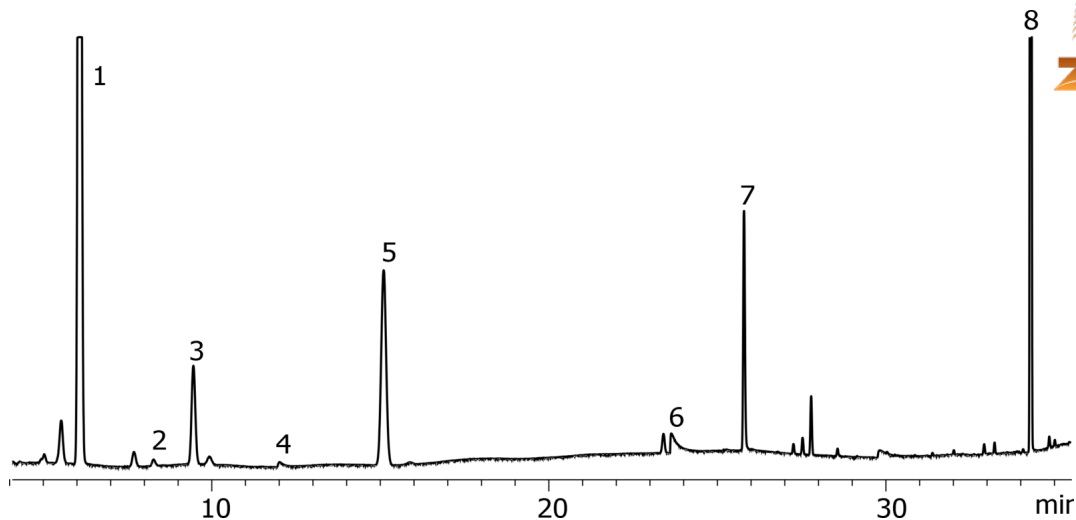