


Cognacs: Hennessey® on ZB-WAXplus**Column:** Zebron™ ZB-WAXplus™, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G013-11**Oven Profile:****Carrier Gas:** Constant Flow Hydrogen, 1 mL/min**Injection:** Split 25.0000000000:1 0 µL @ 210°C**Detection:** Refractive Index (230°C)**Analyst Note:** Static headspace injection of 1ml Hennessey in 20 ml headspace vial at 80 °C for 20 min.**Products used in this application:**

Cognacs: Hennessey[®] on ZB-WAXplus

ANALYTES:

- 1 Acetaldehyde
- 2 Isobutanal
- 3 Ethyl formate
- 4 Ethyl acetate
- 5 Acetal
- 6 Methanol
- 7 Ethanol
- 8 1-Propanol
- 9 Isobutanol
- 10 Methyl-2-butanol
- 11 Methyl-3-butanol
- 12 Ethyl octanoate

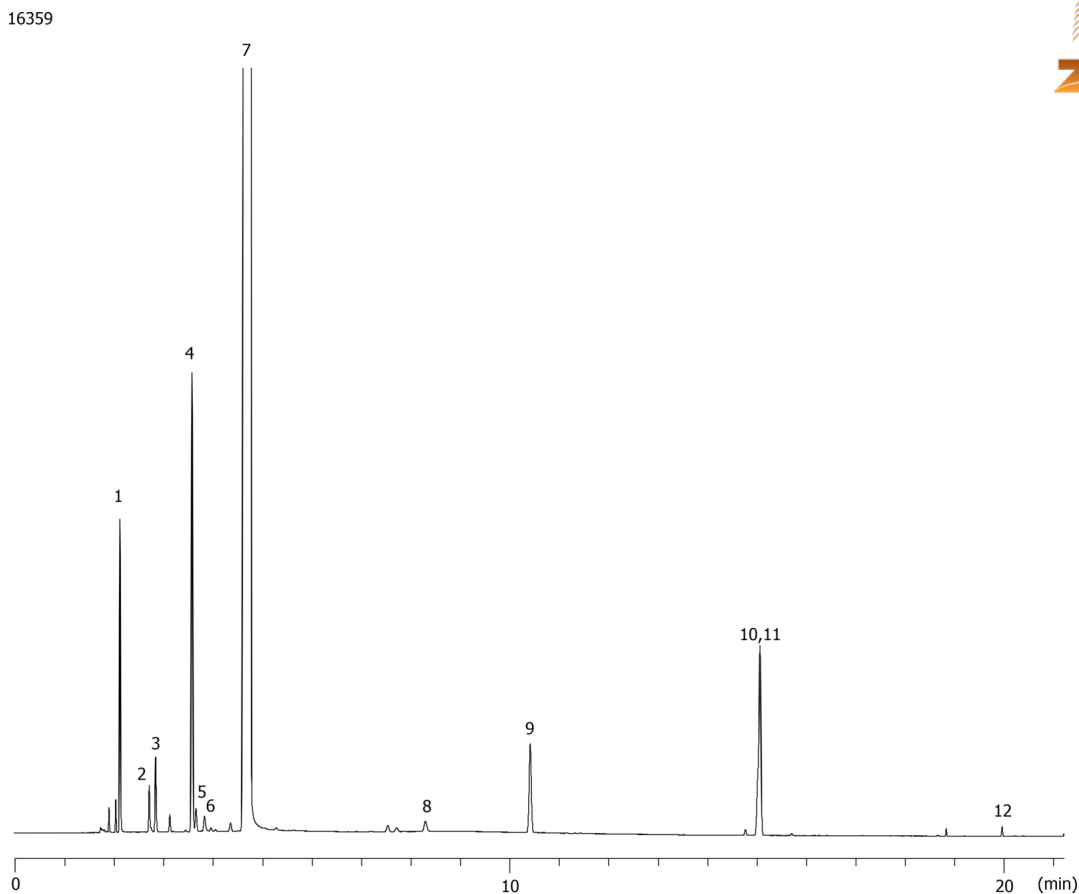