


EPA 8081A Additional Compounds on ZB-5MSi**Column:** Zebron™ ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G018-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 20.000000000:1 0 µL @ 250°C**Detection:** Refractive Index (180°C)**Analyst Note:** Analytes 50ppm in acetone**Products used in this application:**

EPA 8081A Additional Compounds on ZB-5MSi

ANALYTES:

- 1 1,2-Dibromo-3-chloropropane
- 2 Hexachlorocyclopentadiene
- 3 Etridiazole
- 4 Chloroneb
- 5 Propachlor
- 6 Di-allate
- 7 Hexachlorobenzene
- 8 Dichlone
- 9 Chlorothalonil
- 10 Alachlor
- 11 Chlorthal
- 12 gamma-Chlordane
- 13 alpha-Chlordane
- 14 trans-Nonachlor
- 15 Nitrofen
- 16 Chlorobenzilate
- 17 Chloropropylate
- 18 Kepone
- 19 Captafol
- 20 Permethrin
- 21 Permethrin

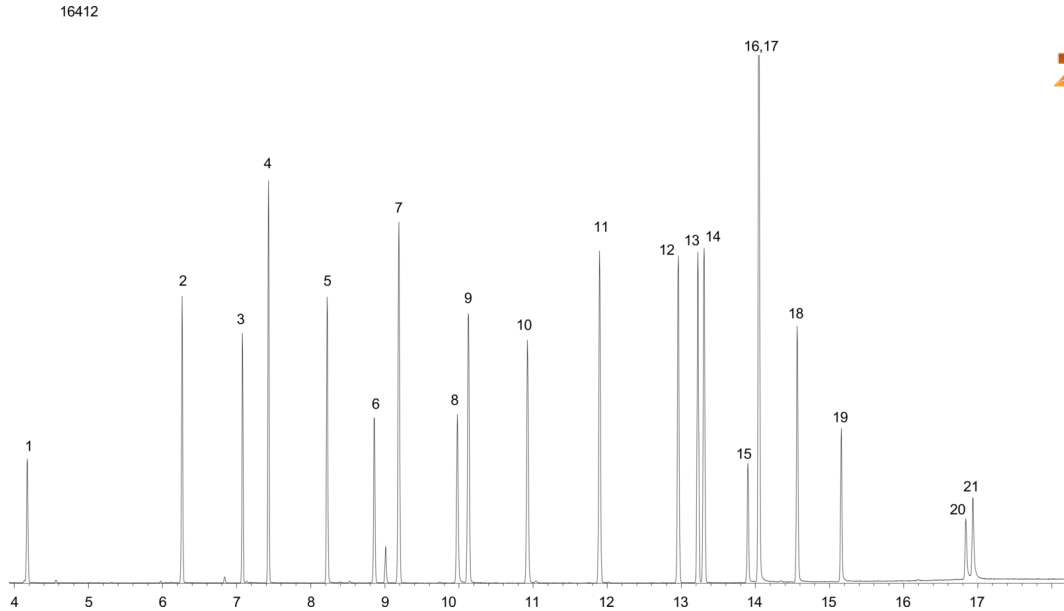