


EPA Method 8082: Aroclor 1232 on ZB-5MSi**Column:** Zebron[™] ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G018-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 70.0000000000:1 0 μ L @ 210°C**Detection:** Refractive Index (310°C)
16416**Products used in this application:****ANALYTES:**

- 1 Tetrachloro-m-xylene
- 2 Decachlorobiphenyl

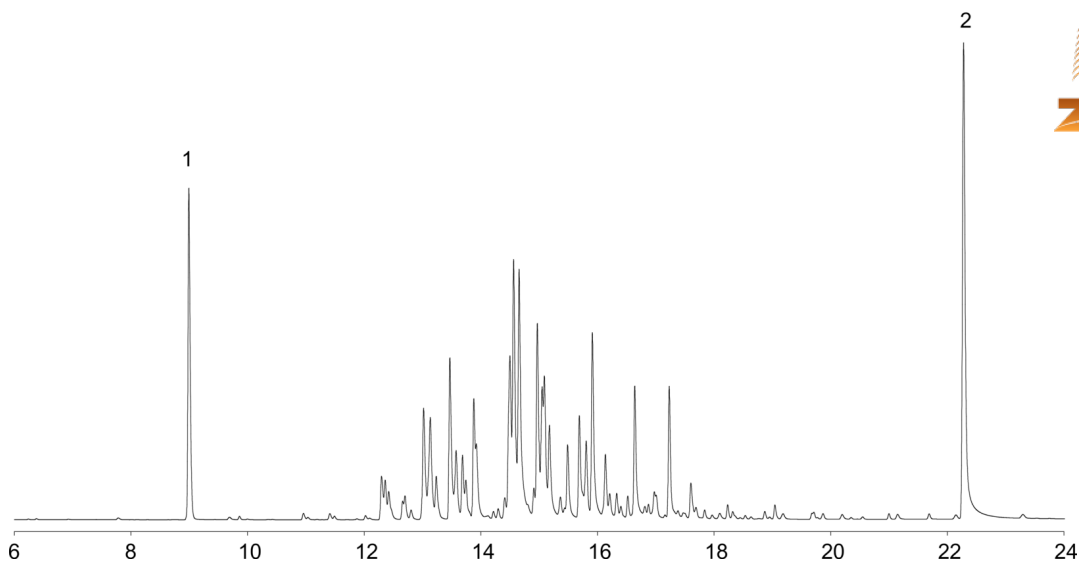