

Acyclovir and Ganciclovir Analysis using Synergi Fusion-RP with MS Detection

Column: Synergi™ 4 µm Fusion-RP 80 Å, LC Column 100 x 2 mm, Ea

Dimensions: 100 x 2 mm ID

Order No: 00D-4424-B0

Elution Type: Isocratic

Eluent A: 10mM Ammonium Formate, pH 3.5

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	97	3	

Flow Rate: 0.4 mL/min


Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 2.0mm, 10/Pk Part No.: AJ0-7556

- Holder Part No.: KJ0-4282

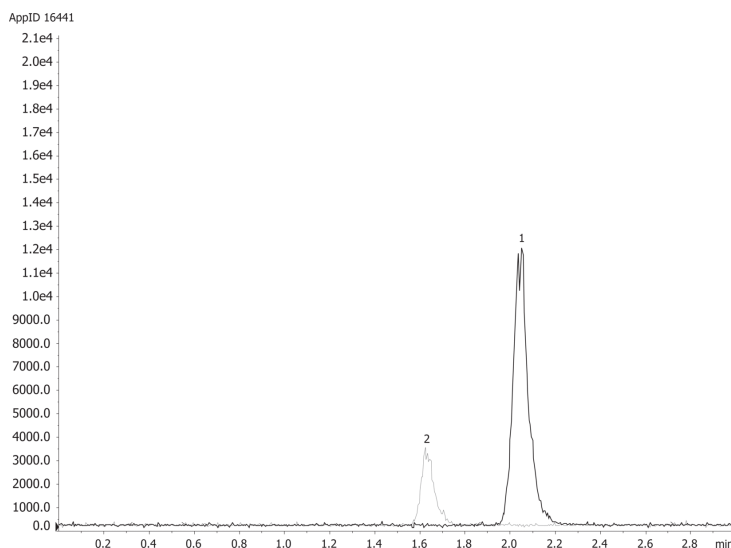

ANALYTES:

1 Acyclovir

Retention Time: 1.62 min

2 Ganciclovir

Retention Time: 2.05 min

