

Improved MS Sensitivity for Bamethan on Luna HILIC

Column: Luna® 3 µm HILIC 200 Å, LC Column 100 x 2 mm, Ea

Dimensions: 100 x 2 mm ID

Order No: 00D-4449-B0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: 100 mM Ammonium Formate, pH 3.2

Gradient	Step No.	Time (min)	Pct A	Pct B
----------	----------	------------	-------	-------

Profile:	1	0	90	10
-----------------	----------	----------	-----------	-----------

Flow Rate: 0.4 mL/min


Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, HILIC 4 x 2.0 mm, 10/Pk Part No.: AJ0-8328

- Holder Part No.: KJ0-4282

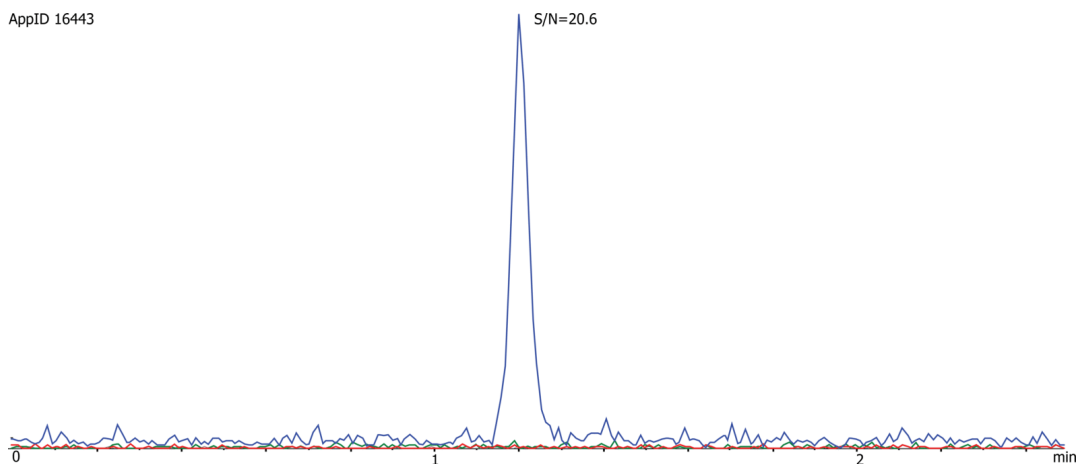

Products used in this application:

**ANALYTES:**

1 Bamethan