


Benzene Derivatives on ZB-5MSi**Column:** Zebron[™] ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G018-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.7 mL/min**Injection:** Splitless 0.000000000:1 0 µL @ 240°C**Detection:** Refractive Index (180°C)**Analyst Note:** 50ppm in methylene chloride**Products used in this application:**

Benzene Derivatives on ZB-5MSi

ANALYTES:

- 1 Methyl methane sulfonate
- 2 Ethyl methane sulfonate
- 3 Benzyl alcohol
- 4 Acetophenone
- 5 Nitrobenzene
- 6 Isophorone
- 7 Safrole
- 8 Isosafrole
- 9 Isosafrole
- 10 1,4-Naphthoquinone
- 11 Dimethyl phthalate
- 12 2,6-Dinitrotoluene
- 13 Dibenzofuran
- 14 2,4-Dinitrotoluene
- 15 Diethyl phthalate
- 16 Azobenzene
- 17 1,3,5-Trinitrobenzene
- 18 Phenacetin
- 19 Dibutyl phthalate
- 20 4-Nitroquinoline-1-oxide
- 21 Butyl benzyl phthalate
- 22 bis(2-Ethylhexyl)phthalate
- 23 Di-n-octyl phthalate

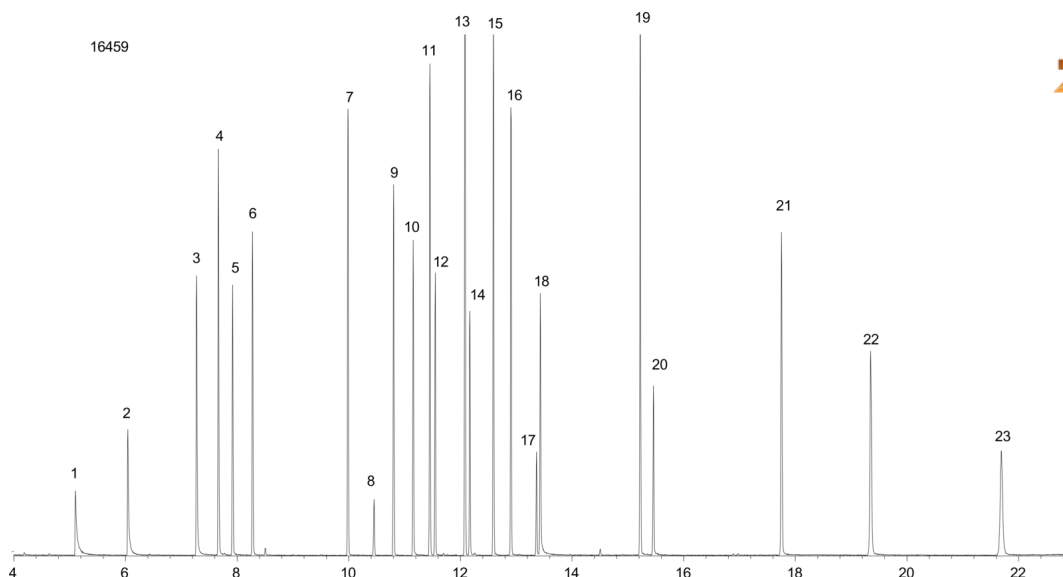