


Drug Screen on ZB-XLB**Column:** Zebron[™] ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.50 μ m, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.5 μ m**Order No:** 7HG-G019-17**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Split 5.0000000000:1 0 μ L @ 250°C**Detection:** Refractive Index (180°C)**Products used in this application:**

Drug Screen on ZB-XLB

ANALYTES:

- 1 Nicotine
- 2 Methylecgonine
- 3 Acetaminophen
- 4 Ibuprofen
- 5 Allobarbitol
- 6 Aprobital
- 7 Butabital
- 8 Phenacetin
- 9 Amobarbitol
- 10 Pentobarbitol
- 11 Meprobamate
- 12 Secobarbitol
- 13 Methyl benzilate
- 14 Benzphetamine
- 15 Dimenhydrinate
- 16 Hexobarbitol
- 17 Caffeine
- 18 Doxylamine
- 19 Phenobarbitol
- 20 Methapyrilene
- 21 Chlorpheniramine
- 22 Brompheniramine
- 23 Cocaine
- 24 Chlorcyclizine
- 25 Codeine
- 26 Diazepam
- 27 Morphine
- 28 Hydrocodone

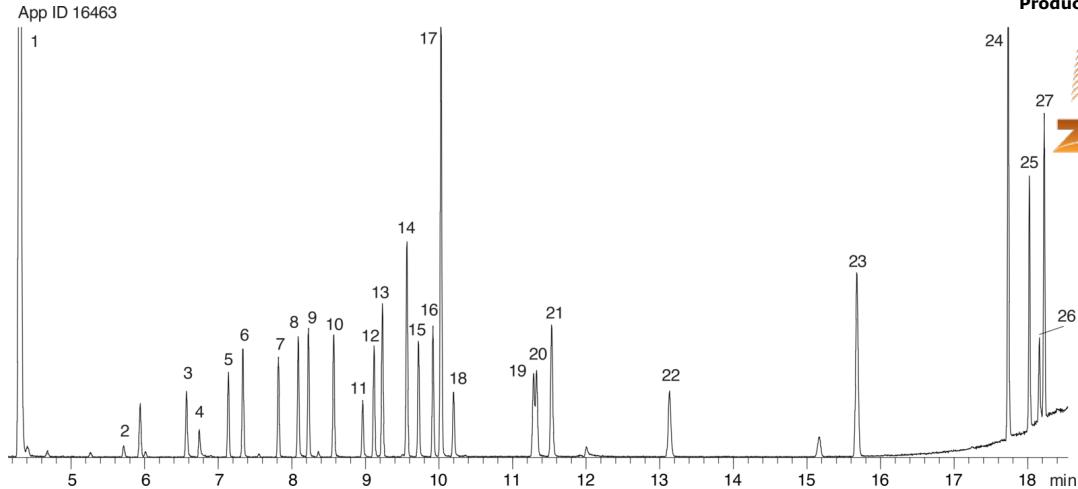