


Polynuclear Aromatic Hydrocarbons (PAHs), EU Mix on ZB-XLB**Column:** Zebtron[™] ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G019-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 5.000000000:1 0 μ L @ 250°C**Detection:** Refractive Index (180°C)**Analyst Note:** 50ppm EU PAH mix in dichloromethane**Products used in this application:**

Polynuclear Aromatic Hydrocarbons (PAHs), EU Mix on ZB-XLB

ANALYTES:

- 1 Benz[a]anthracene
- 2 Cyclopenta(c,d)pyrene
- 3 Chrysene
- 4 5-Methylchrysene
- 5 Benzo[j]fluoranthene
- 6 Benzo[b]fluoranthene
- 7 Benzo[k]fluoranthene
- 8 Benzo[a]pyrene
- 9 Indeno[1,2,3-cd]pyrene
- 10 Dibenz[a,h]anthracene
- 11 Benzo[g,h,i]perylene
- 12 Dibenzo[a,l]pyrene
- 13 Dibenzo[a,e]pyrene
- 14 Dibenzo[a,i]pyrene
- 15 Dibenzo[a,h]pyrene

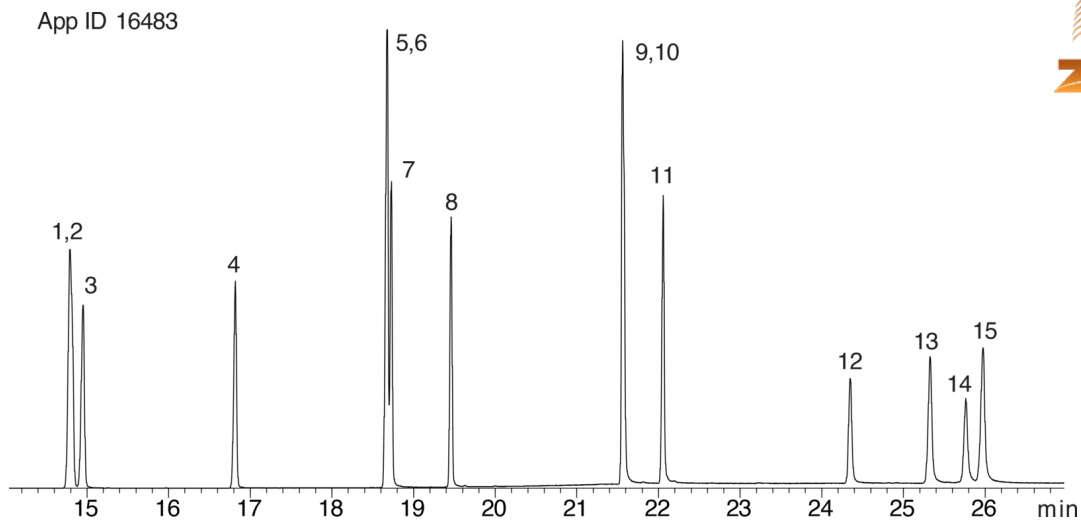