


Polynuclear Aromatic Hydrocarbons (PAHs), EPA 8270 mix on ZB-XLB**Column:** Zebron[™] ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G019-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 5.000000000:1 0 µL @ 250°C**Detection:** Refractive Index (230°C)**Analyst Note:** 50ppm in dichloromethane**Products used in this application:**

Polynuclear Aromatic Hydrocarbons (PAHs), EPA 8270 mix on ZB-XLB

ANALYTES:

- 1 Naphthalene
- 2 2-Methylnaphthalene
- 3 1-Methylnaphthalene
- 4 Acenaphthylene
- 5 Acenaphthene
- 6 Fluorene
- 7 Phenanthrene
- 8 Anthracene
- 9 Fluoranthene
- 10 Pyrene
- 11 Benz[a]anthracene
- 12 Chrysene
- 13 Benzo[b]fluoranthene
- 14 Benzo[k]fluoranthene
- 15 Benzo[a]pyrene
- 16 3-Methylcholanthrene
- 17 Indeno[1,2,3-cd]pyrene
- 18 Dibenzo[a,h]anthracene
- 19 Benzo[g,h,i]perylene

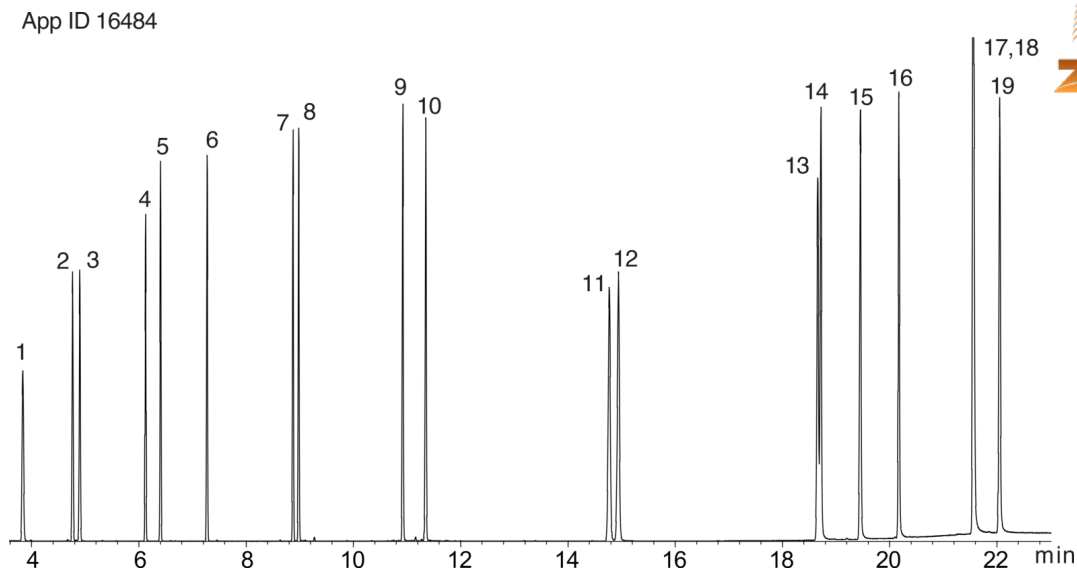