


Glycols in Toothpaste Sample on ZB-WAXplus**Column:** Zebron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G013-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Split 20.0000000000:1 0 μ L @ 225°C**Detection:** Refractive Index (250°C)

Analyst Note: 0.10 mg toothpaste in 1.0 mL of 50-50 ACN/Water. Mixed then centrifuged. 100 μ L supernatant diluted to 400 μ L 50-50 ACN/Water. To inject, transfer 250 μ L to an insert in sample vial with 25 μ L 5mg/mL internal standard. Used a 5 m retention gap pre-column. Glycol standards are in application 16507.
Note - no glycols were found in this toothpaste sample.

**Products used in this application:**

16508

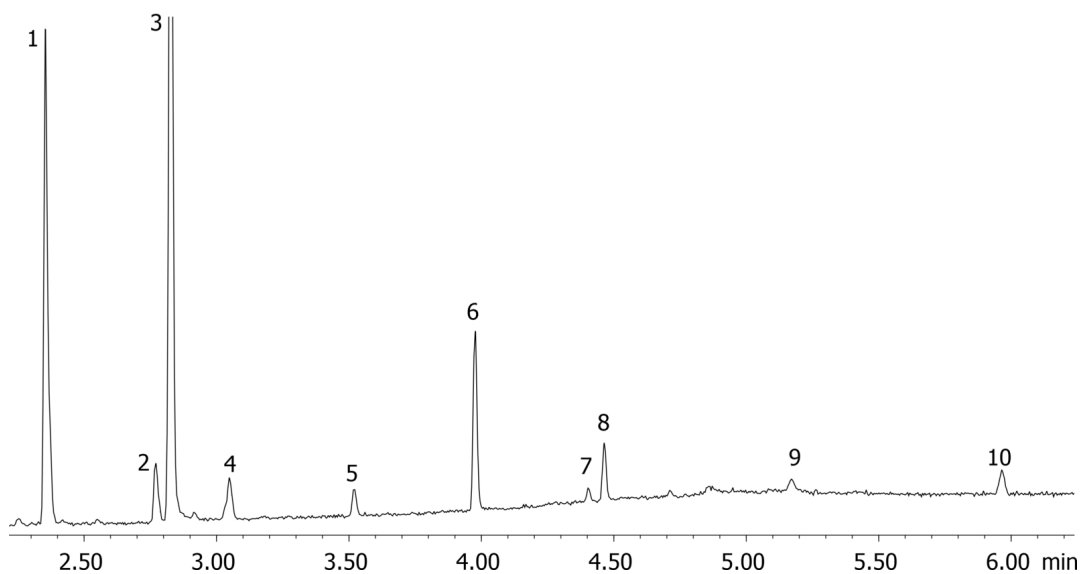

Glycols in Toothpaste Sample on ZB-WAXplus

ANALYTES:

- 1 Levomenthol
- 2 Carvone
- 3 1,3-Propanediol
- 4 Anethole
- 5 Cyclodecane
- 6 3-Phenyl-2-propenal
- 7 Eicosene
- 8 Eugenol
- 9 Glycerin
- 10 o-Methoxycinnamaldehyde

