


Sterols from Lard on ZB-5MS**Column:** Zebron™ ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2.5 mL/min**Injection:** Splitless 0.000000000:1.0 µL @ 320°C**Detection:** Refractive Index (320°C)**Analyst Note:** Sample prep- Saponification, Solid Phase Extraction (SPE), and Derivatization.

16514

**Products used in this application:****ANALYTES:**

- 1 Cholesterol
- 2 Betulin

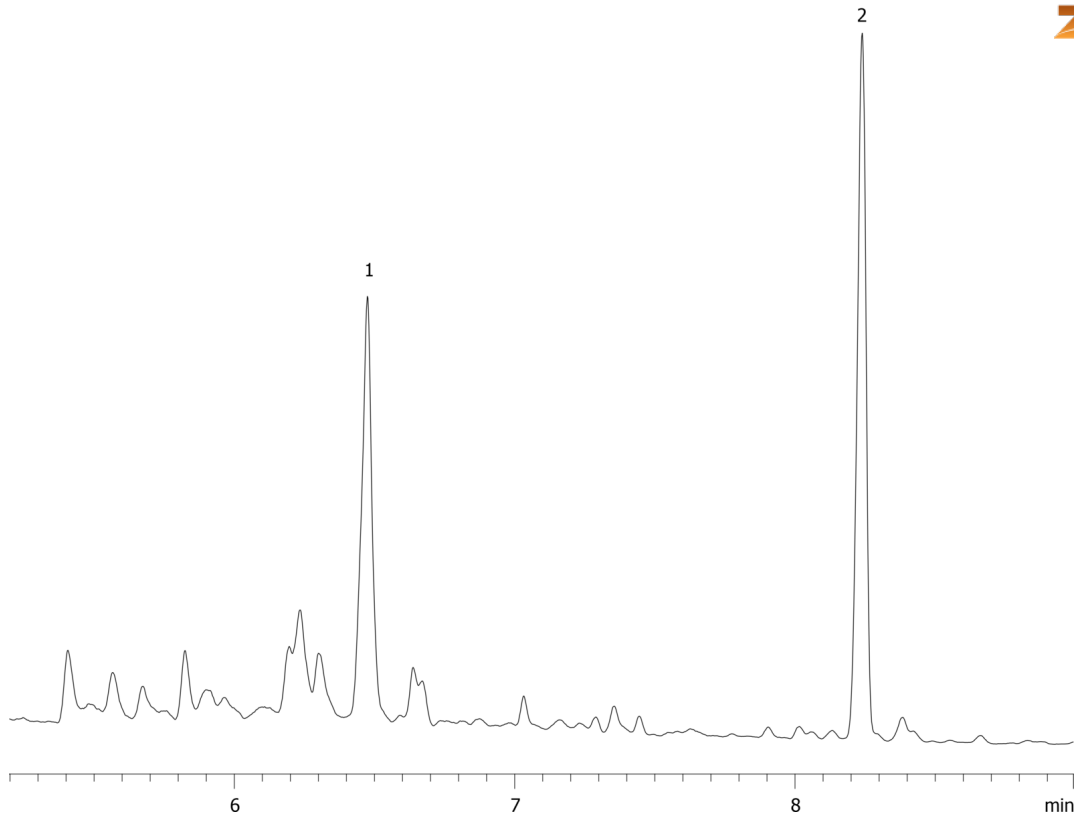