

Organoamines

Column: Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea

Phase: Proprietary Pesticides Phase

Dimensions: 30 meters x 0.25 mm x 0.2 µm

Order No: 7HG-G017-10

Oven Profile:

Carrier Gas: Constant Flow Helium, 1.5 mL/min

Injection: Splitless 0.000000000:1.0 µL @ 240°C

Detection: Refractive Index (320°C)


Analyst Note: Analytes are 50ppm in Dichloromethane



Products used in this application:



16519

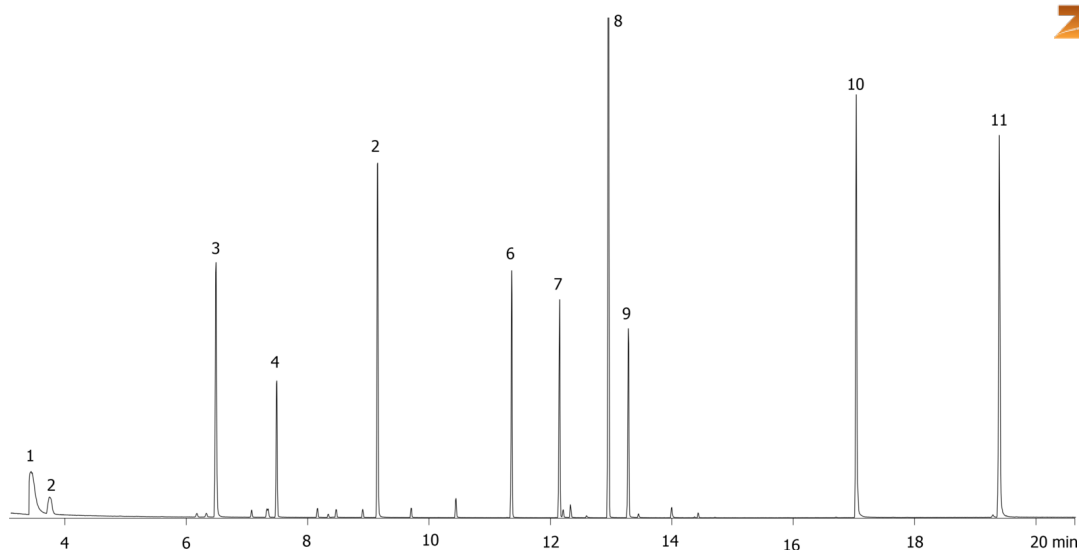

Organoamines

ANALYTES:

- 1 Pyridine
- 2 N-Nitrosodimethylamine
- 3 Aniline
- 4 N-Nitroso-di-n-propylamine
- 5 4-Chloroaniline
- 6 2-Nitroaniline
- 7 3-Nitroaniline
- 8 Diphenylamine
- 9 4-Nitroaniline
- 10 Benzidine
- 11 3,3'-Dichlorobenzidine

