


Organoamines on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 5.000000000:1 0 μ L @ 240°C**Detection:** Refractive Index (230°C)**Analyst Note:** Analytes are 25ppm in methylene chloride**Products used in this application:**

16523

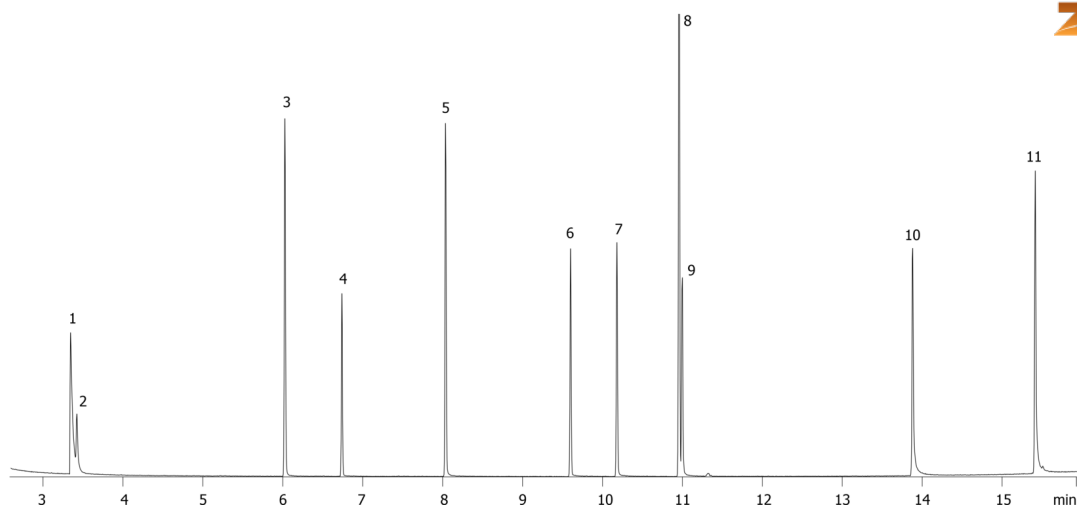

Organoamines on ZB-MR1

ANALYTES:

- 1 Pyridine
- 2 N-Nitrosodimethylamine
- 3 Aniline
- 4 N-Nitrosodi-n-propylamine
- 5 4-Chloroaniline
- 6 2-Nitroaniline
- 7 3-Nitroaniline
- 8 Diphenylamine
- 9 4-Nitroaniline
- 10 Benzidine
- 11 3,3'-Dichlorobenzidine

