


Sterols from Oils and Fats by FID on ZB-MR1**Column:** Zebtron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2.5 mL/min**Injection:** Splitless 0.000000000:1.0 μ L @ 300°C**Detection:** Refractive Index (300°C)**Analyst Note:** Sterols @ 125 ug/ml in 100 μ L pyridine, and 100 μ L BSTFA:TMCS; 99:1, derivatizing agent,**Products used in this application:****ANALYTES:**

- 1 Cholesterol
- 2 Campesterol
- 3 Stigmasterol
- 4 β -Sitosterol
- 5 Betulin

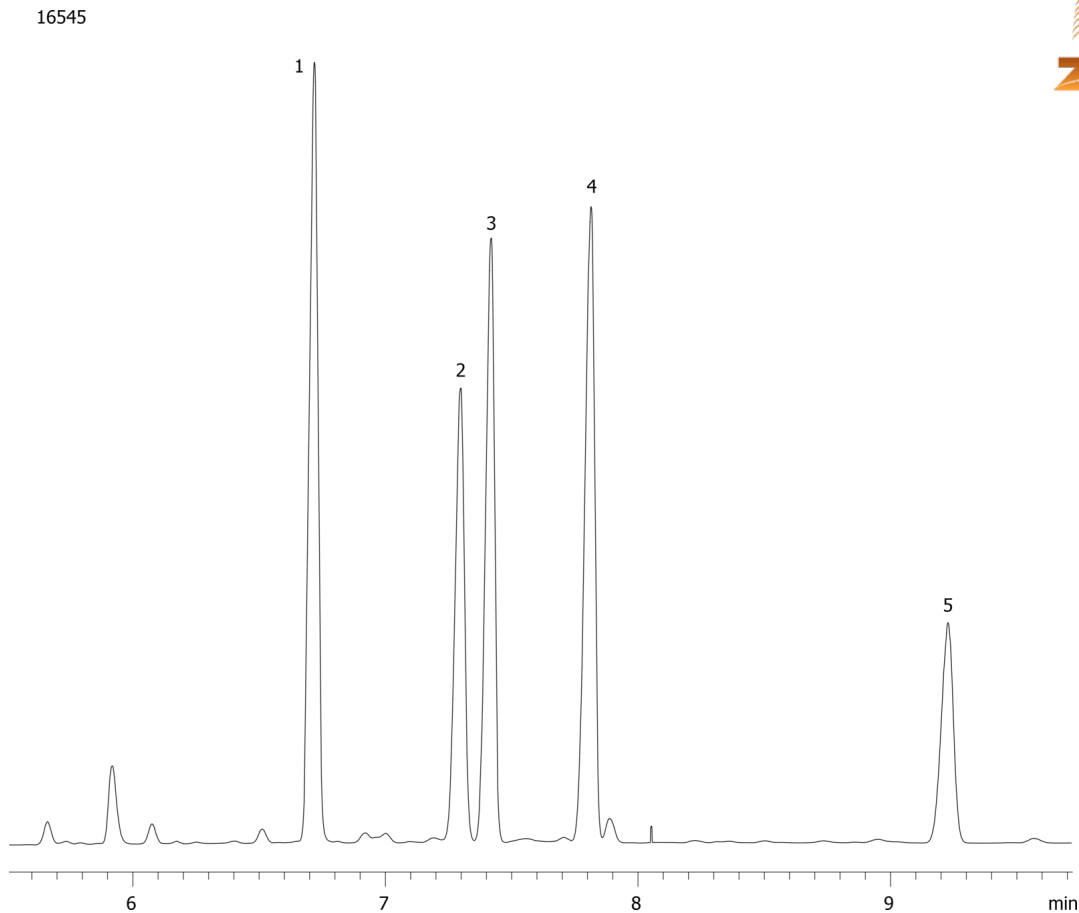