


Organophosphorus pesticides using MSD on ZB-5MSi**Column:** ZebtronTM ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G018-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 20.0000000000:1 0 μ L @ 250°C**Detection:** Refractive Index (180°C)**Analyst Note:** Analytes are 25ppm in isooctane**Products used in this application:**

Organophosphorus pesticides using MSD on ZB-5MSi

ANALYTES:

- 1 Dichlorphos
- 2 Mevinphos
- 3 Demeton
- 4 Ethoprop
- 5 Naled
- 6 Phorate
- 7 Demeton
- 8 Diazinon
- 9 Disulfoton
- 10 Methyl parathion
- 11 Ronnel
- 12 Fenthion
- 13 Chlorpyrifos
- 14 Trichloronate
- 15 Tetrachlorvinphos
- 16 Tokuthion
- 17 Fensulfothion
- 18 Bolstar
- 19 Azinphos-methyl
- 20 Coumaphos

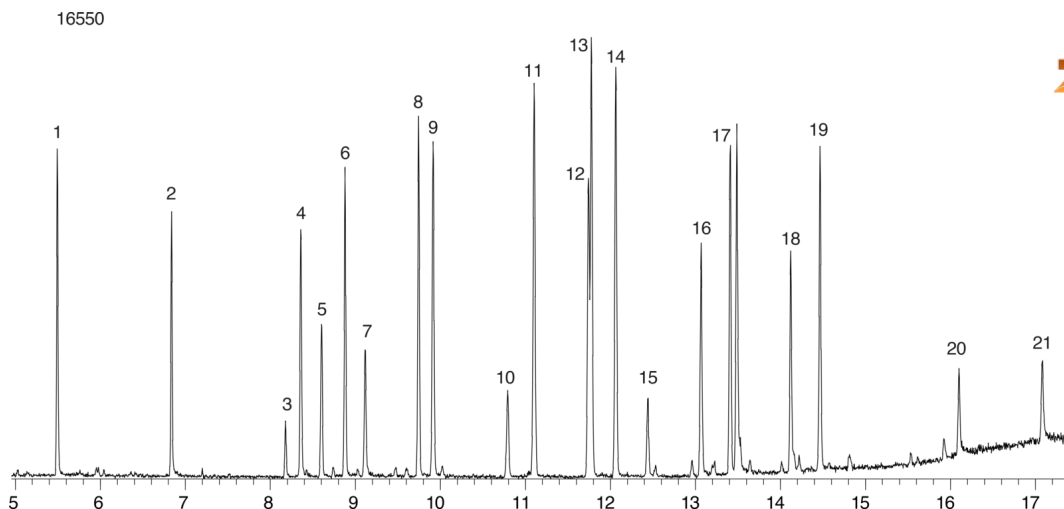