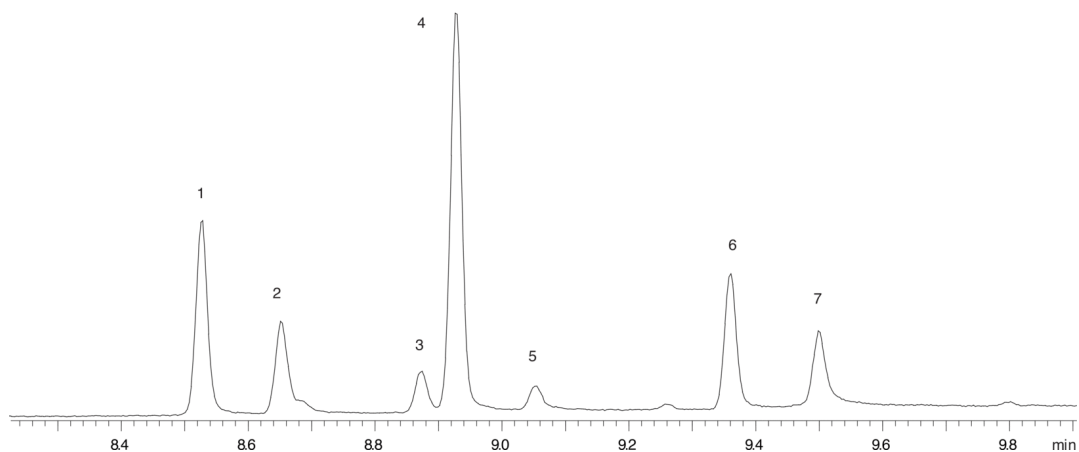


Opiates expanded run by MSD on a ZB-XLB**Column:** Zebtron™ ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G019-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Splitless 5.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (340°C)**Analyst Note:** Analytes are 1ppm in methanol

Sim ions:

(Quant ion*) Codeine (*299,229,162): Norcodeine (*285,215,148): Morphine (*285,215,162): Hydrocodone (*299,242,185): Hydromorphone (*285,229,96): Oxycodone (*315,230,215): Oxymorphone (*301,216,140)

16554

**Products used in this application:****ANALYTES:**

- 1 Codeine
- 2 Norcodeine
- 3 Morphine
- 4 Hydrocodone
- 5 Hydromorphone
- 6 Oxycodone
- 7 Oxymorphone

